



Software Guide



Software Development Kit V3.13

v2.8 Sept 2016

Table of Contents

SECTION 1	INSTALLATION	5
1.1	Technical Support	5
1.2	Rationale	6
1.3	Structure	6
1.4	Key Features	6
1.5	Installation	7
1.5.1	Windows Installation	7
1.5.2	Linux Installation	8
1.6	Getting Started	9
1.6.1	Windows Getting Started	9
1.6.2	Linux Getting Started	11
1.6.3	Microsoft Application Verifier	11
SECTION 2	TUTORIAL	12
2.1	Further Examples	15
2.1.1	Initialize Library and Open Camera	15
2.1.2	Simple Single Frame Acquisition	16
2.1.3	Using a Feature	16
2.1.4	Circular Buffer	17
2.1.7	Pixel Encoding	18
2.1.8	Call-Backs	19
2.1.9	Metadata	20
2.1.10	Binning	22
2.1.11	Acquisition Events	22
2.1.12	Fast Exposure Switching	22
2.2	ESD Recovery	23
SECTION 3	API (APPLICATION PROGRAM INTERFACE)	24
3.1	Overview	24
3.2	Function Listing	24
3.3	API Description	25
3.3.1	Library Initialization	25
3.3.2	Opening a Camera Handle	25
3.3.3	Integer Features	26
3.3.4	Floating Point Features	26
3.3.5	Boolean Features	27
3.3.6	Enumerated Features	27
3.3.7	Command Features	28

3.3.8 String Features.....	29
3.3.9 Buffer Management.....	29
3.3.10 Feature Access Control.....	30
3.3.11 Feature Notifications	31
3.4 Error Codes	33
SECTION 4 FEATURES	36
4.1 Camera Support.....	36
4.2 Feature Reference	36
4.3 Image Format.....	48
4.4 Pixel Encoding.....	50
4.5 Metadata	52
4.6 Area of Interest.....	54
4.7 PixelEncoding and PreAmpGainControl	55
4.8 Sensor Cooling.....	56
4.9 Comparison of SDK2 and SDK3	57
SECTION 5 FUNCTION REFERENCE.....	58
5.1 Function Listing	58
5.1.1 AT_Open	58
5.1.2 AT_Close.....	58
5.1.3 AT_IsImplemented	58
5.1.4 AT_IsReadOnly.....	58
5.1.5 AT_IsWritable.....	59
5.1.6 AT_IsReadable.....	59
5.1.7 AT_RegisterFeatureCallback	59
5.1.8 AT_UnregisterFeatureCallback.....	59
5.1.9 AT_InitialiseLibrary	59
5.1.10 AT_FinaliseLibrary	60
5.1.11 AT_SetInt	60
5.1.12 – AT_GetInt.....	60
5.1.13 – AT_GetIntMax	60
5.1.14 – AT_GetIntMin	60
5.1.15 – AT_SetFloat	60
5.1.16 – AT_GetFloat.....	61
5.1.17 – AT_GetFloatMax	61
5.1.18 – AT_GetFloatMin	61
5.1.19 – AT_SetBool	61
5.1.20 – AT_GetBool.....	62

5.1.21 – AT_Command	62
5.1.22 – AT_SetString	62
5.1.23 – AT_GetString	62
5.1.24 – AT_GetStringMaxLength	62
5.1.25 – AT_SetEnumIndex	63
5.1.26 – AT_SetEnumString	63
5.1.27 – AT_GetEnumIndex	63
5.1.28 – AT_GetEnumCount	63
5.1.29 – AT_GetEnumStringByIndex	63
5.1.30 – AT_IsEnumIndexAvailable	64
5.1.31 – AT_IsEnumIndexImplemented	64
5.1.32 – AT_QueueBuffer	64
5.1.33 – AT_WaitBuffer	64
5.1.34 – AT_Flush	64
SECTION 6 ADDITIONAL LIBRARIES	65
6.1 ATUTILITY	65
6.1.1 – AT_InitialiseUtilityLibrary	65
6.1.2 – AT_FinaliseUtilityLibrary	65
6.1.3 – AT_ConvertBuffer	65
6.1.3.1 Convert Buffer Example	66
6.1.4 – AT_ConvertBufferUsingMetadata	68
6.1.5 – AT_GetWidthFromMetadata	68
6.1.6 – AT_GetHeightFromMetadata	69
6.1.7 – AT_GetStrideFromMetadata	69
6.1.8 – AT_GetPixelEncodingFromMetadata	70
6.1.9 – AT_GetTimeStampFromMetadata	70

SECTION 1 INSTALLATION

1.1 TECHNICAL SUPPORT

If you have any questions regarding the use of this equipment, please contact the representative, from whom your system was purchased, or:

Europe

Andor Technology
7 Millennium Way
Springvale Business Park
Belfast
BT12 7AL
Northern Ireland

Tel. +44 (0) 28 9023 7126
Fax. +44 (0) 28 9031 0792
e-mail: productsupport@andor.com

USA

Andor Technology
425 Sullivan Avenue
Suite # 3
South Windsor
CT 06074
USA

Tel. (860) 290-9211
Fax. (860) 290-9566
e-mail: productsupport@andor.com

Asia-Pacific

Andor Technology (Japan)
7F Ichibancho Central Building
22-1 Ichiban-Cho,
Chiyoda-Ku
Tokyo 102-0082
Japan

Tel. +81 3 3511 0659
Fax. +81 3 3511 0662
e-mail: productsupport@andor.com

China

Andor Technology
Rm 1213
Building B
Luo Ke Time Square
No. 103 Huizhongli
Chaoyang District
Beijing 100101
China

Tel. +86-10-5129-4977
Fax. +86-10-6445-5401
e-mail: productsupport@andor.com

1.2 RATIONALE

Andor SDK Version 3, herein referred to as SDK3, has been designed from the ground up to simplify integration of the Andor camera range into your application. Modern scientific digital cameras have become feature rich devices which can be tailored to the particular application into which they are applied. Andor understands that the integration of the camera is just one component of a larger system solution, and, as such a more consistent and scalable API is required to allow the application developer to both quickly prototype basic acquisition functionality and to provide a clear path to exposing the full feature set.

1.3 STRUCTURE

Install Frame Grabber Card:

- Shut down your PC.
- Install the PCI Express frame grabber card into a free PCI Express slot on your motherboard. Minimum PCIe x8 for Zyla.
- Power on the PC.
- If you are running Windows 7, the drivers will be installed automatically during start up.
- If you are running Windows Vista you may be asked to specify the directory containing the Bitflow drivers, these are located in 'C:\BitFlow SDK 5.60\PlugAndPlay'
- If you are running Windows XP, then the 'Found New Hardware Wizard' dialog will be displayed following the restart.
 - At the 'Can Windows connect to Windows Update to search for software' prompt, select 'No, not this time' and click on the 'Next' button
 - Select 'Install the software automatically' option and click on Install.
 - The Camera Link drivers will now be installed.

1.4 KEY FEATURES

- Simplified API to help reduce development time
- Full access to the current state and limits of camera features
- Support for querying the availability of camera features
- Observer interface to camera features
- Handle parameter in each function to facilitate multiple camera support
- Simple Queue / Wait interface for acquisition buffer management
- Built in software simulated camera (SimCam)

1.5 INSTALLATION

The following sections will describe how to install the software and hardware in order to make your sCMOS camera ready to use.

1.5.1 WINDOWS INSTALLATION

Ensure that you do not install the PCI Express frame grabber card before running any of the software installations.

1. Installation of SDK3 and Camera Link drivers:

[Note: You must have administrator access on your PC to perform the installation.]

- Run the setup.exe file on the cd or from download.
- Select the installation directory or accept the default when prompted by the installer.
- Click on the Install button to confirm and continue with the installation.
- During the installation a number of other windows will pop up as the Camera Link drivers and SDK3 are installed. Click on the Finish button when prompted.

2. Install Frame Grabber Card:

- Shut down your PC.
- Install the PCI Express frame grabber card into a free PCI Express slot on your motherboard.
 - Minimum x4 PCIe for Zyla 3-Tap and Neo 3-Tap
 - Minimum x8 PCIe for Zyla 10-Tap
- Power on the PC.
- **Windows 7:** the drivers are installed automatically during start up.
- **Windows Vista:** you may be asked to specify the directory containing the Bitflow drivers, these are located in 'C:\BitFlow SDK 5.60\PlugAndPlay'
- **Windows XP:** the 'Found New Hardware Wizard' dialog will be displayed following the restart.
 - At the 'Can Windows connect to Windows Update to search for software' prompt, select 'No, not this time' and click on the 'Next' button
 - Select 'Install the software automatically' option and click on Install.
 - The Camera Link drivers will now be installed.

1.5.2 LINUX INSTALLATION

The SDK3 and Camera Link drivers are distributed as a Linux tar file named andor-sdk3-A.B.C.D.tgz where A.B.C.D is the distribution version code.

1. Installation of SDK3 and Camera Link drivers:

- After you have downloaded the tar file you should open a terminal window and change to the download directory.

- Untar the download file by typing:

```
tar -xf andor-sdk3-A.B.C.D.tgz
```

where A.B.C.D is replaced by the version information of the file you have downloaded e.g.

```
tar -xf andor-sdk3-3.1.30005.0.tgz
```

- This creates a sub-folder "andor".
- Change to the "andor" directory and type:

```
sudo ./install_andor
```

- If the install script is unable to determine the platform you will be prompted to enter the platform i.e.

Platform cannot be automatically determined. Please select platform to install:

1. 32-bit

2. 64-bit

3. Exit

Selection:

- Enter appropriate selection, e.g. 2
- The following warning will then be displayed:

This setup will install several libraries into,

/usr/local/lib and

/usr/local/bin

Continue (y/n)?

- Enter "y" to continue.
- If the installation is successful you should see the following:

Bitflow Installation successful

Additional manual configuration required

See the 'BitflowManualConfig.txt' in the 'doc' folder

Andor Installation successful

See the 'doc' directory for further information.

- For CL Cameras: To finish the installation perform the manual steps described in the BitflowManualConfig.txt document. These steps are required to complete the configuration of the cameralink card.
- For USB3 Cameras: To finish the installation perform the manual steps described in the USBManualConfig.txt document.

1.6 GETTING STARTED

This section will demonstrate how to create a basic SDK3 application that will test the software & hardware installation, and will help to confirm communication between PC and camera.

1.6.1 WINDOWS GETTING STARTED

Running the examples that came with the installation

In the installation directory there is an examples directory. In that directory there are two further directories, 'acquisition' and 'serialnumber'. These directories contain the executables and the required DLLs to initialise and communicate with the sCMOS camera.

- The **acquisition** example will initialise and take a single acquisition with the camera and display the counts of the first 20 pixels.
- The **serial number** example will initialise and print out the serial number of the camera.

Creating your own applications

With this installation you can create an application with an Embarcadero or Microsoft compatible compiler. Perform the following steps to create your application.

1. Create a simple console application with either Embarcadero C++ Builder or Microsoft Visual Studio.
2. Add the SDK3 installation directory to the include path for the project. E.g. C:\Program Files\Andor SDK3
3. Add the appropriate library from the SDK3 installation directory to your project.
 - atcore.lib for the Embarcadero compiler
 - atcorem.lib for the Microsoft compiler
4. Copy all the DLL's from the SDK3 installation directory to the directory that the executable is going to run from.
5. Type or copy the code shown below into your projects main file.
6. Compile and run the program. The program should initialise the first camera found and display it's serial number.
7. If the serial number is not displayed then follow the comments in the code listing for hints on tracking down any issues.

NOTE

It is assumed that there is at least 1 sCMOS camera attached. It is acceptable to have no sCMOS cameras attached if using the software simulated camera (SimCam).

```
#include "atcore.h"
#include <iostream>
using namespace std;

int main(int argc, char* argv[])
{
    int i_retCode;
    i_retCode = AT_InitialiseLibrary();
    if (i_retCode != AT_SUCCESS) {
        //error condition, check atdebug.log file
    }
    AT_64 iNumberDevices = 0;
    i_returnCode = AT_GetInt(AT_HANDLE_SYSTEM, L"DeviceCount", &iNumberDevices);
    if (iNumberDevices <= 0) {
        // No cameras found, check all redistributable binaries
        // have been copied to the executable directory or are in the system path
        // and check atdebug.log file
    }
    else {
```

```
AT_H Hndl;
i_retCode = AT_Open(0, &Hndl);
if (i_retCode != AT_SUCCESS) {
    //error condition - check atdebug.log
}

AT_WC szValue[64];
i_retCode= AT_GetString(Hndl, L"SerialNumber", szValue, 64);
if (i_retCode == AT_SUCCESS) {
    //The serial number of the camera is szValue
    wcout << L"The serial number is " << szValue << endl;
}
else {
    //Serial Number feature was not found, check the error code for information
}
AT_Close(Hndl);
}
AT_FinaliseLibrary();

return 0;
}
```

1.6.2 LINUX GETTING STARTED

Running the examples that came with the installation

In the installation directory there is an examples directory. In that directory there are two further directories, 'listdevices' and 'image'. These directories contain source and makefiles. To build the examples change into the appropriate directory and type "make".

1. To run the listdevices example type `./listdevices` from the "listdevices" directory. You should see the following output:
*Found 3 Devices.
Device 0 : DC-152Q-FI
Device 1 : SIMCAM CMOS
Device 2 : SIMCAM CMOS
Press any key and enter to exit.*
2. To run the image example type `./image` from the "image" directory. A bitmap file "image.bmp" should be created in the "image" directory containing a single image acquired from the camera.

Creating your own applications

The example directories contain example source code and makefiles that show how to create your own application. Also see windows example code shown in **Section 1.6.1** Windows Getting Started.

1.6.3 MICROSOFT APPLICATION VERIFIER

It has been found that the SDK3 cannot be used with the Microsoft Application Verifier. This is due to how libusb handles disconnected or inactive devices during initialisation which, though it is later handled within libusb, causes Microsoft Application Verifier to throw an exception.

SECTION 2 TUTORIAL

This section explains how to develop a program to use the API to communicate with a sCMOS camera. It points out the critical parts that the program must have and provides a foundation to build more complicated programs in the future.

The first thing that must be done is to add the appropriate library file to the project. During the SDK3 installation a Borland library (atcore.lib) and a Microsoft library (atcorem.lib) are made available. You should also add the include path of the SDK3 installation directory to your project. The atcore.h file from this directory will need to be included in the main project file.

The very first API call must be AT_InitialiseLibrary, and the very last call must be AT_FinaliseLibrary. These functions will prepare the API for use and free resources when no longer needed.

```
AT_InitialiseLibrary( );
```

```
AT_FinaliseLibrary( );
```

Every function will return an error code when called. It is recommended that a user check every return code before moving on to the next statement. If the AT_InitialiseLibrary function call fails there is no point continuing on with the program as every proceeding function call will also fail.

Every function call will return an Integer and each return code that could possibly be returned is listed in the atcore.h file and documented in **Section 3.4 Error Codes**.

So now the program becomes:

```
int i_returnCode;  
i_returnCode = AT_InitialiseLibrary( );  
if (i_returnCode == AT_SUCCESS) {  
    //continue with program  
}  
i_returnCode = AT_FinaliseLibrary( );  
if (i_returnCode != AT_SUCCESS) {  
    //Error FinaliseLibrary
```

Error Checking

It is highly recommended that you check return codes from every function in case of error. For the purpose of this tutorial the error checking will be kept to a minimum to reduce the length of the program.

The next stage of the program is to get a handle to the camera that is to be controlled. This is done with the AT_Open function and there is a corresponding AT_Close function to release the camera handle. A camera index is passed into the AT_Open function specifying which camera you wish to open.

SimCam

There are two software simulated cameras provided along with the SDK, which allows a user to develop code and see a program running without a hardware camera being connected. If there are any Neo cameras connected then their indexes will start from 0. If the atdevsimcam.dll is in the executable directory then camera indexes after the Neo indexes will reference these software cameras, known as SimCams. We will assume that there is a atdevsimcam.dll and that there is one Neo camera connected, and that we want to control the hardware camera.

Now the program is:

```
int i_returnCode;
AT_H Hndl;
int i_cameraIndex = 0;
i_returnCode = AT_InitialiseLibrary( );
if (i_returnCode == AT_SUCCESS) {
    i_returnCode = AT_Open ( i_cameraIndex, &Hndl );
    if (i_returnCode == AT_SUCCESS) {
        //continue on with program
        //.....
        i_returnCode = AT_Close ( Hndl );
        if (i_returnCode != AT_SUCCESS) {
            // error closing handle
        }
    }
}
i_returnCode = AT_FinaliseLibrary( );
if (i_returnCode != AT_SUCCESS) {
    //Error FinaliseLibrary
}
```

From here on example code will assume that we have successfully got a handle to our camera, following on from the previous code. Now we will add code to the program to modify and view the “ExposureTime” setting of the camera.

It is recommended that any features that need to set or view that they should have their accessibility checked first. Every feature can be checked for being implemented (AT_IsImplemented), writable (AT_IsWritable), readable (AT_IsReadable) and read only (AT_IsReadOnly)

For this tutorial it is assumed that the ‘Exposure Time’ feature is writable and readable, to reduce code complexity.

```
double d_newExposure = 0.02;
i_returnCode = AT_SetFloat ( Hndl, L"ExposureTime", d_newExposure);
if (i_returnCode == AT_SUCCESS) {
    //it has been set
}
```

NOTE

It is recommended that when using the AT_SetFloat functions that an AT_GetFloat is done afterwards to get the actual value that the camera will use as it may not be exactly what was set.

In order to ensure a low noise level in the images we must enable the sensor cooling mechanism. Cooling is off by default, to activate the mechanism use the Boolean feature SensorCooling. Then set the temperature with the TemperatureControl feature and check the status using TemperatureStatus.

```
double temperature = 0;
AT_SetBool(Hndl, L"SensorCooling", AT_TRUE);

int temperatureCount = 0;
AT_GetEnumCount(Hndl, L"TemperatureControl", &temperatureCount);
AT_SetEnumIndex(Hndl, L"TemperatureControl", temperatureCount-1);

int temperatureStatusIndex = 0;
wchar_t* temperatureStatus[256];
AT_GetEnumIndex(Hndl, L"TemperatureStatus", &temperatureStatusIndex);
AT_GetEnumStringByIndex(Hndl, L"TemperatureStatus", temperatureStatusIndex,
temperatureStatus, 256);

//Wait for temperature to stabalise
while(wcscmp(L"Stabilised",temperatureStatus) != 0) {
    AT_GetEnumIndex(Hndl, L"TemperatureStatus", &temperatureStatusIndex);
    AT_GetEnumStringByIndex(Hndl, L"TemperatureStatus", temperatureStatusIndex,
temperatureStatus, 256);
}
```

When setting AccumulateCount, ensure that FrameCount has been set to any value which is a multiple of the value of AccumulateCount. For example, if AccumulateCount has been set to 4, valid values for FrameCount are 4, 8, 12, 16 etc. This ensures that each accumulation has the correct number of frames; failure to do so will result in a timeout error during the acquisition.

All other settings we will take as default.

This stage shows how to get data from a single acquisition. In preparation the program must allocate memory to store the acquired image in. To get the size of memory to be declared use the integer "ImageSizeBytes" feature. No error checking will be shown in the following example to help with clarity- but it is recommended that in final code all return codes are checked.

Buffer Byte Alignment

The memory declared to hold the acquisition should be aligned on an 8 byte boundary, this helps with system performance and also prevent any alignment fault that could occur. The pucAlignedBuffer variable in the following code shows how to enforce 8 byte alignment.

```
AT_64 ImageSizeBytes;
AT_GetInt( Hndl, L"ImageSizeBytes", &ImageSizeBytes);
//cast so that the value can be used in the AT_QueueBuffer function
int i_imageSize = static_cast<int>(ImageSizeBytes);

unsigned char* uc_Buffer = NULL;
gblp_Buffer = new unsigned char[i_imageSize+8]; // Add 8 to allow data alignment

// Adjust pointer so that it falls on an 8-byte boundary
unsigned char* pucAlignedBuffer = reinterpret_cast<unsigned char*>(
(reinterpret_cast<unsigned long>( gblp_Buffer ) + 7 ) & ~0x7);
```

The next stage is to let SDK know what memory to use for the upcoming acquisition. This is done with the AT_QueueBuffer API function call. Multiple buffers can be queued to the SDK before an acquisition starts if you are acquiring a sequence of frames. For now we will assume that only one frame is required.

```
AT_QueueBuffer(Hndl, pucAlignedBuffer, ImageSizeBytes);
```

Now start the acquisition with the Command "AcquisitionStart". The command to stop the acquisition is "AcquisitionStop".

```
AT_Command(Hndl, L"AcquisitionStart");  
//get the data  
AT_Command(Hndl, L"AcquisitionStop");
```

When the acquisition has been started there is a function `AT_WaitBuffer` that can be used to put the calling thread to block until the current image has been captured and is ready for the program to use. A time out value in milliseconds is also specified to `AT_WaitBuffer` to force it to return if the acquisition has not occurred in that time frame.

```
unsigned char* pBuf;  
int BufSize;  
AT_WaitBuffer(Hndl, &pBuf, &BufSize, 10000);
```

It is vital to check the return code from the `AT_WaitBuffer` function before processing the returned buffer. The `pBuf` return from this function should be the same pointer to the one that was queued in the `AT_QueueBuffer` function. If all has been successful, then the data in the array pointed to by `pBuf` will contain the acquired image. The "Pixel Encoding" feature should then be checked to confirm the format of the datastream. See the **Features** Section for a more complete explanation.

After the acquisition is complete any buffers queued up with the `AT_QueueBuffer` command but not yet returned from `AT_WaitBuffer`, need to be released by the SDK or else the next acquisition will use them.-Failure to call the `AT_Flush` function after an acquisition will not only result in both queues not being cleared, but may lead to undefined behaviour.

```
AT_Flush(Hndl);
```

The final code for this tutorial can be seen in section "Code Listing for Tutorial" in the Appendix.

2.1 FURTHER EXAMPLES

2.1.1 INITIALIZE LIBRARY AND OPEN CAMERA

```
//InitialiseLibrary must be the first function call made by an application before  
accessing other functions  
AT_InitialiseLibrary();  
  
//Declare an Andor Device Handle for referencing device later  
AT_H Handle;  
  
//Open the first device (Device 0)  
AT_Open(0, &Handle);  
  
//Close the device when finished using it  
AT_Close(Handle);  
  
//Call FinaliseLibrary when all access to API is complete  
AT_FinaliseLibrary();
```

2.1.2 SIMPLE SINGLE FRAME ACQUISITION

```
AT_InitialiseLibrary();
AT_H Handle;
AT_Open(0, &Handle);

//Set the exposure time for this camera to 10 milliseconds

AT_SetFloat(Handle, L"ExposureTime", 0.01);

//Get the number of bytes required to store one frame
AT_64 ImageSizeBytes;
AT_GetInt(Handle, L"ImageSizeBytes", &ImageSizeBytes);

int BufferSize = static_cast<int>(ImageSizeBytes);

//Allocate a memory buffer to store one frame
unsigned char* UserBuffer = new unsigned char[BufferSize];

//Pass this buffer to the SDK
AT_QueueBuffer(Handle, UserBuffer, BufferSize);

//Start the Acquisition running
AT_Command(Handle, L"AcquisitionStart");

//Sleep in this thread until data is ready, in this case set
//the timeout to infinite for simplicity
unsigned char* Buffer;
AT_WaitBuffer(Handle, &Buffer, &BufferSize, AT_INFINITE);

//Stop the acquisition
AT_Command(Handle, L"AcquisitionStop");
AT_Flush(Handle);

//Application specific data processing goes here..

//Free the allocated buffer
delete [] UserBuffer;

AT_Close(Handle);
AT_Close(Handle);
AT_FinaliseLibrary();
```

2.1.3 USING A FEATURE

```
AT_InitialiseLibrary();
AT_H Handle;
AT_Open(0, &Handle);

AT_BOOL Implemented;
//To determine if Exposure time is implemented by the camera
AT_IsImplemented(Handle, L"ExposureTime", &Implemented);

AT_BOOL ReadOnly;
//To determine if Exposure time a Read Only Feature
AT_IsReadOnly(Handle, L"ExposureTime", &ReadOnly);

if (Implemented==AT_TRUE) {
```



```
//Get the Limits for Exposure Time
double Min, Max;
AT_GetFloatMin(Handle, L"ExposureTime", &Min);
AT_GetFloatMax(Handle, L"ExposureTime", &Max);

//Get the current accessibility
AT_BOOL Writable, Readable;
AT_IsWritable(Handle, L"ExposureTime", &Writable);
AT_IsReadable(Handle, L"ExposureTime", &Readable);

if (Readable==AT_TRUE) {
    //To get the current value of Exposure time in
    //microseconds
    double ExposureTime;
    AT_GetFloat(Handle, L"ExposureTime", &ExposureTime);
}

if (Writable==AT_TRUE) {
    //To set the value of Exposure Time to 10
    //microseconds
    AT_SetFloat(Handle, L"ExposureTime", 0.00001);
}

}

AT_Close(Handle);
AT_FinaliseLibrary();
```

2.1.4 CIRCULAR BUFFER

```
AT_InitialiseLibrary();
AT_H Handle;
AT_Open(0, &Handle);

AT_SetFloat(Handle, L"ExposureTime", 0.01);

AT_64 ImageSizeBytes;
AT_GetInt(Handle, L"ImageSizeBytes", &ImageSizeBytes);

int BufferSize = static_cast<int>(ImageSizeBytes);

//Declare the number of buffers and the number of frames interested in
int NumberOfBuffers = 10;
int NumberOfFrames = 100;

//Allocate a number of memory buffers to store frames
unsigned char** AcqBuffers = new unsigned char*[NumberOfBuffers];
unsigned char** AlignedBuffers = new unsigned char*[NumberOfBuffers];
for (int i=0; i < NumberOfBuffers; i++) {
    AcqBuffers[i] = new unsigned char[BufferSize + 7];
    AlignedBuffers[i] = reinterpret_cast<unsigned char*>((reinterpret_cast<unsigned
long>(AcqBuffers[i% NumberOfBuffers]) + 7) & ~7);
}

//Pass these buffers to the SDK
for(int i=0; i < NumberOfBuffers; i++) {
    AT_QueueBuffer(Handle, AlignedBuffers[i], BufferSize);
}
```

```
//Set the camera to continuously acquires frames
AT_SetEnumString(Handle, L"CycleMode", L"Continuous");

//Start the Acquisition running
AT_Command(Handle, L"AcquisitionStart");

//Sleep in this thread until data is ready, in this case set
//the timeout to infinite for simplicity
unsigned char* pBuf;
int BufSize;
for (int i=0; i < NumberOfFrames; i++) {
    AT_WaitBuffer(Handle, &pBuf, &BufSize, AT_INFINITE);

    //Application specific data processing goes here..

    //Re-queue the buffers
    AT_QueueBuffer(Handle, AlignedBuffers[i*NumberOfBuffers], BufferSize);
}

//Stop the acquisition
AT_Command(Handle, L"AcquisitionStop");
AT_Flush(Handle);

//Application specific data processing goes here..

//Free the allocated buffer
for (int i=0; i < NumberOfBuffers; i++) {
    delete[] AcqBuffers[i];
}
delete[] AlignedBuffers;
delete[] AcqBuffers;

AT_Close(Handle);
AT_FinaliseLibrary();
```

2.1.7 PIXEL ENCODING

```
#define EXTRACTLOWPACKED(SourcePtr) ( (SourcePtr[0] << 4) + (SourcePtr[1] & 0xF) )
#define EXTRACTHIGHPACKED(SourcePtr) ( (SourcePtr[2] << 4) + (SourcePtr[1] >> 4) )

AT_InitialiseLibrary();
AT_H Handle;
AT_Open(0, &Handle);

//Set the pixel Encoding to the desired settings Mono12Packed Data
AT_SetEnumString(Handle, L"PixelEncoding", L"Mono12Packed");

AT_64 ImageSizeBytes;
AT_GetInt(Handle, L"ImageSizeBytes", &ImageSizeBytes);

int BufferSize = static_cast<int>(ImageSizeBytes);
unsigned char* UserBuffer = new unsigned char[BufferSize];

AT_QueueBuffer(Handle, UserBuffer, BufferSize);
AT_Command(Handle, L"AcquisitionStart");

unsigned char* Buffer;
AT_WaitBuffer(Handle, &Buffer, &BufferSize, AT_INFINITE);
```

```
//Stop the acquisition
AT_Command(Handle, L"AcquisitionStop");
AT_Flush(Handle);
//Unpack the 12 bit packed data
for (int i=0; i < BufferSize; i+=3) {
    AT_64 LowPixel = EXTRACTLOWPACKED(Buffer);
    AT_64 HighPixel = EXTRACTHIGHPACKED(Buffer);

    //Application specific data processing goes here..

    Buffer += 3;
}

delete [] UserBuffer;

AT_Close(Handle);
AT_FinaliseLibrary();
```

2.1.8 CALL-BACKS

//This example will demonstrate how to setup, register and unregister a call-back
 //Tests of the correct call-back context and a count of the number of updates received
 are provided

```
AT_H Handle;
int g_iCallbackCount = 0;
int g_iCallbackContext = 0;

int AT_EXP_CONV Callback(AT_H Hndl, const AT_WC* Feature, void* Context)
{
    //Application specific call-back handling should go here
    g_iCallbackCount++;
    g_iCallbackContext = *reinterpret_cast<int*>(Context);
    return AT_CALLBACK_SUCCESS;
}

int main(int argc, char* argv[])
{
    AT_InitialiseLibrary();
    AT_Open(0, &Handle);

    //Set the call-back context, context values can be defined on per application basis
    int i_callbackContext = 5;

    //Reset the call-back count
    //Only required for the purposes of this example to show the call-back has been
    received
    g_iCallbackCount = 0;

    //Register a call-back for the given feature
    AT_RegisterFeatureCallback(Handle, L"PixelReadoutRate", Callback,
    (void*)&i_callbackContext);

    //Set the feature in order to trigger the call-back
    AT_SetEnumIndex(Handle, L"PixelReadoutRate", 0);

    // Application specific code should go here
    //For this example we shall check that the call-back has been successful
```

```

if (g_iCallbackCount==0 || g_iCallbackContext != i_callbackContext) {
    //Deal with failed call-back
}

//Unregister the call-back, no more updates will be received
AT_UnregisterFeatureCallback(Handle, L"PixelReadoutRate", Callback,
(void*)&i_callbackContext);

AT_Close(Handle);
AT_FinaliseLibrary();
}

```

2.1.9 METADATA

The following examples illustrate enabling metadata when taking acquisitions, however, they do not show how to correctly parse the results. For more details on parsing metadata refer to the `metadataparser` example located in the examples folder of the SDK3 installation directory. Alternatively functions are available in `ATUtility` to help you retrieve metadata. See **Section 6.1 ATUTILITY**.

Toggle Metadata

```

...
int i_returnCode;
AT_BOOL i_metadataEnabled;
i_returnCode = AT_GetBool(Hndl, L"MetadataEnable", &i_metadataEnabled);
if (i_returnCode != AT_SUCCESS) {
    // Problem getting state of metadata
}

i_metadataEnabled = ( i_metadataEnabled == AT_FALSE ) ? AT_TRUE : AT_FALSE;

i_returnCode = AT_SetBool(Hndl, L"MetadataEnable", i_metadataEnabled);
if (i_returnCode != AT_SUCCESS) {
    // Problem enabling metadata
}
...

```

Parsing Timestamp

```

...
#define LENGTH_FIELD_SIZE 4
#define CID_FIELD_SIZE 4
#define TIMESTAMP_FIELD_SIZE 8

// Setup environment. (See example in appendix).
...
// Take an acquisition. (See example in appendix).
unsigned char* puc_image;
...

// Get image size.
AT_64 ImageSizeBytes;
AT_GetInt(Hndl, L"ImageSizeBytes", &ImageSizeBytes);
int i_imageSize = static_cast<int>(ImageSizeBytes);

// Move to end of image. This is assuming reading metadata right to left.
unsigned char* puc_metadata = puc_image + i_imageSize;

// Extract length field from end.
int i_length = *(reinterpret_cast<int*>(puc_metadata - LENGTH_FIELD_SIZE));

```

```
// Move to start of timestamp information. Note that this example assumes that
// timestamps is the first metadata block. Real code should search for the
// timestamp CID
int i_offset = LENGTH_FIELD_SIZE + CID_FIELD_SIZE + TIMESTAMP_FIELD_SIZE;
AT_64 au64_timestamp = *(reinterpret_cast<andoru64*>(puc_metadata - i_offset));

// Tear down environment. (See example in appendix).
...
```

Parsing Frame Info

```
...
#define LENGTH_FIELD_SIZE 4
#define CID_FIELD_SIZE 4
#define FRAMEINFO_FIELD_SIZE 8

// Setup environment. (See example in appendix).
...
// Take an acquisition. (See example in appendix).
unsigned char* puc_image;
...

// Get image size.
AT_64 ImageSizeBytes;
AT_GetInt(Hndl, L"ImageSizeBytes", &ImageSizeBytes);
int i_imageSize = static_cast<int>(ImageSizeBytes);

// Move to end of image. This is assuming reading metadata right to left.
unsigned char* puc_metadata = puc_image + i_imageSize;

// Extract length field from end.
int i_length = *(reinterpret_cast<int*>(puc_metadata - LENGTH_FIELD_SIZE));

// Move to start of frame info information. Note that this example assumes that
// frame info is the first metadata block. Real code should search for the
// frame info CID
int i_offset = LENGTH_FIELD_SIZE + CID_FIELD_SIZE + FRAMEINFO_FIELD_SIZE;
AT_64 au64_frameInfo = *(reinterpret_cast<andoru64*>(puc_metadata - i_offset));

AT_64 au64_stride = static_cast<AT_64>(au64_frameInfo & 0xFFFF);
au64_frameInfo = au64_frameInfo >> 16;
AT_64 au64_pixelEncoding = static_cast<AT_64>(au64_frameInfo & 0xFF);
au64_frameInfo = au64_frameInfo >> 16;
AT_64 au64_aoiWidth = static_cast<AT_64>(au64_frameInfo & 0xFFFF);
au64_frameInfo = au64_frameInfo >> 16;
AT_64 au64_aoiHeight = static_cast<AT_64>(au64_frameInfo & 0xFFFF);

// Tear down environment. (See example in appendix).
```

2.1.10 BINNING

Configure Binning for SCMOS

```
int i_returnCode;
i_returnCode = AT_SetEnumString(Hndl, L"AOIBinning", L"3x3");
if (i_returnCode != AT_SUCCESS) {
    // Problem setting binning
}
```

2.1.11 ACQUISITION EVENTS

```
// Setup environment. (See example in appendix).
...

// Enable Exposure End Events
int returnCode = AT_SetEnumString(Handle, L"EventSelector",
                                   L"ExposureEndEvent");
if (AT_SUCCESS != returnCode) {
    // Problem selecting Exposure End Event
}
returnCode = AT_SetBool(Handle, L"EventEnable", AT_TRUE);
if (AT_SUCCESS != returnCode) {
    // Problem enabling Exposure End Event
}

// Register a Callback for Exposure End Events. (See Callback example)
AT_RegisterFeatureCallback(Handle, L"ExposureEndEvent", Callback, NULL);

// Take an acquisition. (See example in appendix).
...

// Disable Exposure End Event
returnCode = AT_SetBool(Handle, L"EventEnable", AT_FALSE);
if (AT_SUCCESS != returnCode) {
    // Problem disabling Exposure End Event
}

AT_UnregisterFeatureCallback(Handle, L"ExposureEndEvent", Callback, NULL);

// Tear down environment. (See example in appendix).
```

2.1.12 FAST EXPOSURE SWITCHING

The exposure time can be changed during an acquisition without first stopping the camera to apply the new time. However, some limitations apply to the fast-switchable exposure range. Thus, some exposure times cannot be set whilst an acquisition is running, and a restart will be required. To determine if a restart is required, establish the min and max values using the following commands:

```
AT_GetFloatMin(handle, L"ExposureTime", &minExposureTime);
AT_GetFloatMax(handle, L"ExposureTime", &maxExposureTime);
```

If the exposure time to be set is within the min and max boundaries, no restart will be required, and can be set during the acquisition. Values outside this range cannot be set whilst the acquisition is running.

2.2 ESD RECOVERY

An effective solution and recovery from an ESD strike to the ZYLA USB3 camera relies on the user application's ability to respond to SDK error codes and callbacks correctly.

The following table outlines the corrective actions that must be taken by the user application to recover from an ESD strike that has led to loss of frames and errors as a result.

The solution consists of 3 parts:

1. Zyla FX3 microcode version 1.09.0 or greater
2. SDK3 upgrade version 3.12.30000 or greater
3. Higher-level software application recommendations

Refer to table below for details of scenarios where user software application action is recommended.

Acquisition Type	Possible ESD Strike Effect (SDK3 application)	Recommended Detection	Recommended Action
None (idle)	None, camera will be automatically re-acquired and remain in the idle state.	Register a callback on the CameraPresent feature	None
Software Triggered	The software trigger may be lost or frame data may be corrupted or dropped. AT_WaitBuffer will timeout for one or more of the frames that have been dropped.	Register a callback on the CameraPresent feature	When camera has finished restarting (CameraPresent is true), abort and restart the acquisition. N.B. it is not sufficient to send another software trigger.
Externally (TTL) Triggered	Frame data may be corrupted or dropped. AT_WaitBuffer will timeout for one or more of the frames that have been dropped.	Register a callback on the CameraPresent feature	When camera has finished restarting (CameraPresent is true), abort and restart the acquisition.
Internally Triggered (free-running)	Frame data may be corrupted or dropped but may not time out. A gap will be noticeable in the camera's metadata (hardware timestamps).	Register a callback on the CameraPresent feature	When camera has finished restarting (CameraPresent is true), abort and restart the acquisition.

SECTION 3 API (APPLICATION PROGRAM INTERFACE)

3.1 OVERVIEW

The SDK3 API can be divided into several sets of functions, each controlling a particular aspect of camera control. There are sections in the API for opening a handle to a camera, for buffer management and for accessing the features that every camera exposes. Each feature that a camera exposes to the user has a particular type that represents how that feature is controlled. The feature types are:

- Integer
- Floating Point
- Boolean
- Enumerated
- Command
- String

For example:

- Exposure Time feature: Floating Point
- Acquisition Start feature: Command

Each of these feature types, the management of multiple cameras and buffer management are described in the sections below. The character type used by the API is a 16 bit wide character defined by the AT_WC type, which is used to represent all feature names, enumerated options and string feature values.

Wide Characters

An example of converting wide character strings to char strings can be found in the appendix.

3.2 FUNCTION LISTING

```
int AT_InitialiseLibrary();
int AT_FinaliseLibrary();

int AT_Open(int DeviceIndex, AT_H* Handle);
int AT_Close(AT_H Hndl);

typedef int (*FeatureCallback)(AT_H Hndl, AT_WC* Feature, void* Context);
int AT_RegisterFeatureCallback(AT_H Hndl, AT_WC* Feature, FeatureCallback EvCallback,
void* Context);
int AT_UnregisterFeatureCallback(AT_H Hndl, AT_WC* Feature, FeatureCallback EvCallback,
void* Context);

int AT_IsImplemented(AT_H Hndl, AT_WC* Feature, AT_BOOL* Implemented);
int AT_IsReadOnly(AT_H Hndl, AT_WC* Feature, AT_BOOL* ReadOnly);
int AT_IsReadable(AT_H Hndl, AT_WC* Feature, AT_BOOL* Readable);
int AT_IsWritable(AT_H Hndl, AT_WC* Feature, AT_BOOL* Writable);

int AT_SetInt(AT_H Hndl, AT_WC* Feature, AT_64 Value);
int AT_GetInt(AT_H Hndl, AT_WC* Feature, AT_64* Value);
int AT_GetIntMax(AT_H Hndl, AT_WC* Feature, AT_64* MaxValue);
int AT_GetIntMin(AT_H Hndl, AT_WC* Feature, AT_64* MinValue);

int AT_SetFloat(AT_H Hndl, AT_WC* Feature, double Value);
int AT_GetFloat(AT_H Hndl, AT_WC* Feature, double* Value);
int AT_GetFloatMax(AT_H Hndl, AT_WC* Feature, double* MaxValue);
int AT_GetFloatMin(AT_H Hndl, AT_WC* Feature, double* MinValue);
```



```
int AT_SetBool(AT_H Hndl, AT_WC* Feature, AT_BOOL Value);
int AT_GetBool(AT_H Hndl, AT_WC* Feature, AT_BOOL* Value);

int AT_SetEnumIndex(AT_H Hndl, AT_WC* Feature, int Value);
int AT_SetEnumString(AT_H Hndl, AT_WC* Feature, AT_WC* String);
int AT_GetEnumIndex(AT_H Hndl, AT_WC* Feature, int* Value);
int AT_GetEnumCount(AT_H Hndl, AT_WC* Feature, int* Count);
int AT_IsEnumIndexAvailable(AT_H Hndl, AT_WC* Feature, int Index, AT_BOOL* Available);
int AT_IsEnumIndexImplemented(AT_H Hndl, AT_WC* Feature, int Index, AT_BOOL* Implemented);
int AT_GetEnumStringByIndex(AT_H Hndl, AT_WC* Feature, int Index, AT_WC* String, int StringLength);

int AT_Command(AT_H Hndl, AT_WC* Feature);

int AT_SetString(AT_H Hndl, AT_WC* Feature, AT_WC* Value);
int AT_GetString(AT_H Hndl, AT_WC* Feature, AT_WC* Value, int StringLength);
int AT_GetStringMaxLength(AT_H Hndl, AT_WC* Feature, int* MaxStringLength);

int AT_QueueBuffer(AT_H Hndl, AT_U8* Ptr, int PtrSize);
int AT_WaitBuffer(AT_H Hndl, AT_U8** Ptr, int* PtrSize, unsigned int Timeout);
int AT_Flush(AT_H Hndl);
```

3.3 API DESCRIPTION

3.3.1 LIBRARY INITIALIZATION

The first API function call made by any application using SDK3 must be:

```
AT_InitialiseLibrary()
```

This allows SDK to setup its internal data structures and to detect any cameras that are attached. AT_InitialiseLibrary takes no parameters.

Before your application closes or when you no longer wish to access the API you should call the function:

```
AT_FinaliseLibrary()
```

This cleans up any data structures held internally by SDK.

3.3.2 OPENING A CAMERA HANDLE

To access the features provided by a camera and to acquire images you must first open a handle. A camera handle, represented by the data type AT_H, is a reference to the particular camera that you wish to control and is passed as the first parameter to most other functions in the SDK. In multi-camera environments the handle becomes particularly useful as it allows cameras to be controlled simultaneously in a thread safe manner. To open a handle to a camera you should pass the index of the camera that you wish to access, to the function:

```
AT_Open(int DeviceIndex, AT_H* Handle)
```

The handle will be returned in the Handle parameter which is passed by address. To open the first camera you should pass a value of 0 to the DeviceIndex parameter, for the second camera pass a value of 1 etc.

Once you have finished with the camera it should be closed using the function:

```
AT_Close(AT_H Handle)
```

The only parameter to this function is the handle of the camera that you wish to release.

System Handle

There are some features of the system that are not connected to a specific camera but are global properties. For example, the Device Count feature stores a count of the number of devices that are currently connected. To access these features you do not need to open a handle to a camera, instead you should use the system handle represented by the constant `AT_HANDLE_SYSTEM`. You do not need to retrieve this handle using the `AT_Open` function; it is predefined and can be used immediately after the `AT_InitialiseLibrary` function has completed.

Thread Safety

SDK3 is thread safe when accessing different cameras on different threads. However you should not use multiple threads at once to access the features of the same camera. The exception to this is that the `AT_WaitBuffer` and `AT_QueueBuffer` functions can be called simultaneously on separate threads, with the same camera handle, during acquisition.

3.3.3 INTEGER FEATURES

Integer features are those that can be represented by a single integer value. For example, the number of images in the sequence that you wish to acquire is represented by the Integer feature `FrameCount`. To set an Integer feature use the function:

```
AT_SetInt(AT_H Hndl, AT_WC* Feature, AT_64 Value)
```

The first parameter is the handle to the camera that is exposing the desired feature; the second parameter is a wide character string indicating the name of the feature that you wish to modify. The full list of feature strings is available in the Feature Reference section. The third parameter contains the value that you want to assign to the feature. The function will return a value indicating whether the function successfully applied the value. The **Error Codes (Section 3.4 Error Codes)** lists the possible error codes that can be returned from the Integer Type functions.

To get the current value for an Integer feature use the function:

```
AT_GetInt(AT_H Hndl, AT_WC* Feature, AT_64 * Value)
```

The first two parameters are the camera handle and the feature name, the same as those passed to the `AT_SetInt` function. The third parameter is the address of the variable into which you want to store the Integer value.

Integer features can sometimes be restricted in the range of values that they can be set to. This range of possible values can be determined by using the functions:

```
AT_GetIntMax(AT_H Hndl, AT_WC* Feature, AT_64 * MaxValue)
AT_GetIntMin(AT_H Hndl, AT_WC* Feature, AT_64 * MinValue)
```

These functions work similarly to the `AT_GetInt` function except that the third parameter returns either the highest allowable value or the lowest allowable value. Using these maximum and minimum values you can check which values are allowed by `AT_SetInt` without having to monitor its return code. This can be useful, for example, when you wish to limit the range of possible values that the user can enter in a GUI application. Note that the maximum and minimum of an Integer feature may change if other dependent features values are modified, for example, the maximum frame rate will decrease as the exposure time is increased. You can use the feature notification mechanism described in a later section to find out when this happens.

3.3.4 FLOATING POINT FEATURES

Floating Point type features work in a similar way to Integer features, in that they have a Set function and Get function and GetMin and GetMax functions. Floating Point features represent those features that are expressed with a

value that contains a decimal point. As an example, the Exposure Time feature is exposed through the Floating Point functions. The functions are:

```
AT_SetFloat(AT_H Hndl, AT_WC* Feature, double Value)
AT_GetFloat(AT_H Hndl, AT_WC* Feature, double * Value)
AT_GetFloatMax(AT_H Hndl, AT_WC* Feature, double * MaxValue)
AT_GetFloatMin(AT_H Hndl, AT_WC* Feature, double * MinValue)
```

The first parameter to each of these functions is the camera handle, the second parameter is the name of the feature and the third parameter contains either the value that you wish to set or the address of a variable that will return the current value, maximum or minimum of the feature. The list of possible error codes is described in the **Error Codes Section 3.4 Error Codes**. Note that the maximum and minimum of a Floating Point feature may change if other dependent features values are modified, you can use the feature notification mechanism described in later section to find out when this happens.

3.3.5 BOOLEAN FEATURES

Boolean features can only be set to one of two possible values, representing the logical states true and false. True is represented by the value AT_TRUE and false by the value AT_FALSE. An example of a boolean feature is the Sensor Cooling feature which can be used to switch the cooler on the camera on or off. To change the state of a boolean feature use the function:

```
AT_SetBool(AT_H Hndl, AT_WC* Feature, AT_BOOL Value)
```

The first parameter is a handle to the camera being used, the second parameter is the string descriptor of the feature to change and the third parameter is the value. So to enable a boolean feature pass a value of AT_TRUE, to disable a boolean feature, pass a value of AT_FALSE in the third parameter. To retrieve the current state of a boolean feature, use the function:

```
AT_GetBool(AT_H Hndl, AT_WC* Feature, AT_BOOL* Value)
```

For this function the third parameter contains the address of the variable into which you want the state stored. A value of AT_FALSE means the feature is disabled, while a value of AT_TRUE means the feature is enabled.

3.3.6 ENUMERATED FEATURES

Enumerated features are used to represent those features that can be assigned one value from a set of possible options. For example, the triggering mode that you wish to use with the camera is set using the TriggerMode enumerated feature. The triggering mode setting can be chosen from a number of options, for example, internal, external or external start. The enumerated feature functions allow you to:

- Determine how many options are available
- Select which option you wish to use
- Retrieve a human readable representation of each option.

Enumerated options can be set either by their text value or by index, using the functions:

```
AT_SetEnumIndex(AT_H Hndl, AT_WC* Feature, int Index)
AT_SetEnumString(AT_H Hndl, AT_WC* Feature, AT_WC* String)
```

The first function changes the current item to the one that lies at the position specified by the Index parameter, where an Index of 0 is the first item. The second function lets you specify the string descriptor for the particular option that you wish to use. String Descriptors for all features are described in the Feature Reference **Section 4.2 Feature Reference**. As for all feature access functions the first two parameters are the camera handle and the string descriptor of the feature that you wish to modify. The choice of which function to use will depend on your particular application and they can both be used in the same program.

Enumerated Indexes

The particular index that maps to an enumerated option may be different across SDK versions and across different cameras. To ensure best compatibility for your application you should use strings wherever possible and avoid assuming that a specific option is found at a particular index.

To find out which option is currently selected for an enumerated feature, you can use the function:

```
AT_GetEnumIndex(AT_H Hndl, AT_WC* Feature, int* Value)
```

The third parameter is the address of the variable where you want the currently selected index stored.

To find out how many options there are available for the feature, use the function:

```
AT_GetEnumCount(AT_H Hndl, AT_WC* Feature, int* Count)
```

The third parameter, on return, will contain the number of possible options. If you attempt to select an option using `AT_SetEnumIndex` with an index either below zero or above or equal to this count an error will be returned.

You can retrieve the string descriptor for any option by calling the function:

```
AT_GetEnumStringByIndex(AT_H Hndl, AT_WC* Feature, int Index, AT_WC* String,
                        int StringLength)
```

The third parameter is the index of the option that you want to receive the descriptor for, the fourth parameter is a user allocated buffer to receive the descriptor and the fifth parameter is the length of the allocated buffer.

Enumerated Index Availability

In some situations one or more of the options listed for an enumerated feature may be either permanently or temporarily unavailable. An option may be permanently unavailable if the camera does not support this option, or temporarily unavailable if the current value of other features do not allow this option to be selected. To find out which options are available you can use the functions:

```
AT_IsEnumIndexAvailable(AT_H Hndl, AT_WC* Feature, int Index, AT_BOOL* Available)
AT_IsEnumIndexImplemented(AT_H Hndl, AT_WC* Feature, int Index,
                          AT_BOOL* Implemented)
```

Both functions take the index of the option that you want to interrogate in the third parameter. The first function determines if the option is only temporarily unavailable, the second function determines if the feature is permanently unavailable. The fourth parameter returns the availability status, a value of `AT_FALSE` means unavailable, and a value of `AT_TRUE` means the option is available. If you try to select an option that is unavailable using either of the set functions then an error will be returned.

3.3.7 COMMAND FEATURES

Command features are those that represent a single action. For example, to start the camera acquiring you will use the Command feature Acquisition Start. These commands do not require any extra parameters and are simply called by passing the string descriptor of the command to the function:

```
AT_Command(AT_H Hndl, AT_WC* Feature)
```

The function call is blocking so when the function returns, the action is complete.

3.3.8 STRING FEATURES

String features are those that can be represented by 1 or more characters. An example of a String feature is the Serial Number of the camera. In many cases these features are read only but if they are writable, you can set the value using the function:

```
AT_SetString(AT_H Hndl, AT_WC* Feature, AT_WC* Value)
```

The first parameter contains the camera handle, the second parameter is the string descriptor of the feature, and the third parameter is the character string that you want to assign to the feature.

To retrieve the value of a String feature use the function:

```
AT_GetString(AT_H Hndl, AT_WC* Feature, AT_WC* Value, int StringLength)
```

In this case the third parameter should be a caller allocated character string that will be used to receive the string. The fourth parameter is the length of the caller allocated buffer. To determine what length of string that should be allocated to receive the string value you can use the function:

```
AT_GetStringMaxLength(AT_H Hndl, AT_WC* Feature, int* MaxStringLength)
```

The maximum length of the String feature will be returned in the third parameter.

3.3.9 BUFFER MANAGEMENT

SDK maintains two queues for each camera, which are used to manage the transfer of image data to the application. Both queues operate in a First-in-First-out (FIFO) basis and are used to store the addresses of blocks of memory allocated by the application.

Queuing

The first queue, the input queue, which is written into by the application and read from by the SDK, is used to store the memory buffers that have not yet been filled with image data. This queue is accessed using the function:

```
AT_QueueBuffer(AT_H Hndl, AT_U8 * Ptr, int PtrSize)
```

The first parameter contains the handle to the camera. The second parameter is the address of an application allocated buffer large enough to store a single image. The PtrSize parameter should contain the size of the buffer that is being queued. The required size of the buffer can be obtained by reading the value of the ImageSizeBytes integer feature. The AT_QueueBuffer function can be called multiple times with different buffers, to allow a backlog of buffers to be stored by the SDK. By doing this the SDK can be copying image data into these buffers while the application is processing previous images.

Buffer Alignment

Any buffers queued to the SDK using the AT_QueueBuffer function should have their address aligned to an 8-byte boundary. The examples shown later in this manual demonstrate how this can be done if your compiler does not do this for you automatically when creating the buffer.

Waiting

The second queue is the output queue and is used to store the application defined buffers after they have had images copied into them. In this case the SDK adds buffers to this queue which can then be retrieved by the application. This application can retrieve the processed buffers from this queue by using the function:

```
AT_WaitBuffer(AT_H Hndl, AT_U8 ** Ptr, int* PtrSize, unsigned int Timeout)
```

The AT_WaitBuffer function will retrieve the next buffer from the output queue and return the address in the second parameter; the size of the buffer will also be returned in the PtrSize parameter. As both of these parameters are

outputs from the function, they are passed in by address. If there are no buffers currently in the output queue, the `AT_WaitBuffer` function will put the calling thread to sleep until a buffer arrives. The thread will sleep until either a buffer arrives or the time specified by the `Timeout` parameter expires. The `Timeout` parameter is specified in milliseconds and can be any value between 0 and the constant `AT_INFINITE`. If a value of zero is used then the function will simply test the output queue for available buffers and return immediately. If the value is `AT_INFINITE`, then the function will sleep indefinitely until data arrives at the output queue; any value in between will be used as a millisecond timeout for the function.

Flushing

The input and output queues are not automatically flushed when an acquisition either completes normally or is stopped prematurely. Any buffers remaining in the input queue will be used during the next acquisition and any buffers in the output queue are still available to be retrieved by the application. If you wish to clear the two queues at any time then you should call the function below. If this function is not called after an acquisition is complete, then the remaining buffers will be used the next time an acquisition is started and may lead to undefined behaviour.

```
AT_Flush(AT_H Hndl)
```

Acquisition Control

See the descriptions of the Acquisition Start and Acquisition Stop command features for information on running an acquisition on the camera.

3.3.10 FEATURE ACCESS CONTROL

The individual access rights of features can be determined using a set of functions that apply to all features, independent of their type. The four access characteristics of a feature are:

- Whether a feature is implemented by a camera.
- Whether a feature is read only.
- Whether a feature can currently be read.
- Whether a feature can currently be modified.

The first two access rights are permanent characteristics of the feature, the second two access rights may change during the running of the program. For example, if other features are modified in such a way as to affect this feature. If a feature is not implemented by a camera then any attempt to access that feature will return the error code `AT_ERR_NOTIMPLEMENTED`. Any attempt to modify a read only feature will result in the error code `AT_ERR_READONLY`. If a feature cannot be currently read then any attempt to get the current value will return the `AT_ERR_NOTREADABLE` error code and any attempt to modify a value that cannot currently be written to will return the error code `AT_ERR_NOTWRITABLE`.

To determine if a feature has been implemented use the function:

```
AT_IsImplemented(AT_H Hndl, AT_WC* Feature, AT_BOOL* Implemented)
```

The first two parameters are, as usual, the handle to the camera and the string descriptor of the feature. The third parameter is an output parameter which returns with a value indicating whether the feature is implemented or not. If `Implemented` contains the value `AT_TRUE`, on return from the function then the feature is implemented, if it returns with the value `AT_FALSE`, then the feature is not implemented.

To determine if a feature is read only, use the function:

```
AT_IsReadOnly(AT_H Hndl, AT_WC* Feature, AT_BOOL* ReadOnly)
```


This function works in a similar way to `AT_IsImplemented`, that is, if the `ReadOnly` parameter returns with the value `AT_TRUE` then the feature is read only, a value of `AT_FALSE` indicates that it can be modified. `SerialNumber` is an example of a feature that is read only.

To determine if a feature is currently readable, use the function:

```
AT_IsReadable(AT_H Hndl, AT_WC* Feature, AT_BOOL* Readable)
```

The `Readable` parameter indicates whether the feature is currently readable in the same manner as the `ReadOnly` parameter to the function `AT_IsReadOnly`.

To determine if a feature is currently writable use the function:

```
AT_IsWritable(AT_H Hndl, AT_WC* Feature, AT_BOOL* Writable)
```

An example of the use of this function is to allow a GUI application to disable access to features when they cannot be modified, for example, whilst an acquisition is running.

3.3.11 FEATURE NOTIFICATIONS

Sometimes a feature may change its value or its other characteristics, not as a direct result of the user modifying the feature, but indirectly through modification of a separate feature. For example, if the `Trigger Mode` feature is set to `External trigger`, then the `Frame Rate` feature's writable access characteristic will be disabled (see `Feature Access Control` above), as the frame rate is now controlled by the rate at which the external trigger is applied, and not by the application setting.

To allow the application to receive notification when this type of indirect change occurs, there are functions provided in the API that allow the application to create a callback function and attach it to a feature. Whenever the feature changes in any way, this callback will be triggered, allowing the application to carry out any actions required to respond to the change. For example, if an application provides a GUI interface that allows users to modify features, then the callback can update the GUI with any changes. This facilitates use of the Observer design pattern in your application. The callback will also be triggered if the feature is modified directly by the application.

The definition of the callback function implemented by the application should be in the format:

```
int AT_EXP_CONV MyCallback(AT_H Hndl, const AT_WC* Feature, void* Context)
{
    // Perform action
}
```

There are three parameters sent to the function that allow the application to determine the reason for the call-back. The first parameter indicates which camera caused the call-back. By using this parameter you can make use of the same call-back function for multiple cameras. The second parameter holds the string descriptor of the feature that has been modified in some way, and allows the same call-back function to be used with multiple features. The final parameter is an application defined context parameter that was passed in as a parameter at the time that the call-back function was registered. The `Context` parameter is not parsed in any way by the SDK and can be used to store any information that the application wishes.

Note that the `AT_EXP_CONV` modifier must be present and ensures that the correct calling convention is used by the SDK.

To register the call-back function, use the function:

```
AT_RegisterFeatureCallback(AT_H Hndl, AT_WC* Feature, FeatureCallback EvCallback, void* Context)
```

The first parameter is the camera handle, the second parameter is the string descriptor of the feature that you wish to receive notifications for. The third parameter is your call-back function that you have defined as described above and

the fourth parameter in the Context parameter that will be passed to the call-back each time it is called. Whenever the application no longer requires notifications for a particular feature, it should release the call-back by calling the function:

```
AT_UnregisterFeatureCallback(AT_H Hndl, AT_WC* Feature, FeatureCallback EvCallback,  
void* Context)
```

The same parameters should be passed to this function as were passed to the AT_RegisterFeatureCallback.

You need to register a call-back individually for each feature that you are interested in, but the same call-back function can be used for all or some features, or, a separate call-back function can be provided for each feature.

Notes on implementing call-backs

A call-back should complete any work required in the minimal amount of time as it holds up the thread that caused the call-back. If possible the application should delegate any work to a separate application thread if the action will take a significant amount of time.

The call-back function should not attempt to modify the value of any feature as this can cause lockup.

3.4 ERROR CODES

Find below the available return codes and their values for each feature type and the buffer control functions.

Device Connection	Description
AT_SUCCESS (0)	Function call has been successful
AT_ERR_NONINITIALISED (1)	Function called with an uninitialized handle
AT_ERR_CONNECTION (10)	Error connecting to or disconnecting from hardware
AT_ERR_INVALIDHANDLE (12)	Invalid device handle passed to function
AT_ERR_NULL_HANDLE (21)	Null device handle passed to function
AT_ERR_NOMEMORY (37)	No memory has been allocated for the current action
AT_ERR_DEVICEINUSE (38)	Function failed to connect to a device because it is already being used

String Feature	Description
AT_SUCCESS (0)	Function call has been successful
AT_ERR_NOTIMPLEMENTED (2)	Feature has not been implemented for the chosen camera
AT_ERR_READONLY (3)	Feature is read only
AT_ERR_NOTWRITABLE (5)	Feature is currently not writable
AT_ERR_NOTREADABLE (4)	Feature is currently not readable
AT_ERR_EXCEEDEDMAXSTRINGLENGTH (9)	String value provided exceeds the maximum allowed length
AT_ERR_NULL_FEATURE (20)	NULL feature name passed to function
AT_ERR_NULL_READABLE_VAR (23)	Readable not set
AT_ERR_NULL_WRITABLE_VAR (25)	Writable not set
AT_ERR_NULL_ISAVAILABLE_VAR (31)	Available not set
AT_ERR_NULL_VALUE (28)	NULL value returned from function
AT_ERR_NULL_STRING (29)	NULL string returned from function
AT_ERR_NULL_MAXSTRINGLENGTH (32)	Max string length is NULL
AT_ERR_INVALIDHANDLE (12)	Invalid device handle passed to function
AT_ERR_NOMEMORY (37)	No memory has been allocated for the current action
AT_ERR_COMM (17)	An error has occurred while communicating with hardware

Integer Feature	Description
AT_SUCCESS (0)	Function call has been successful
AT_ERR_OUTOFRANGE (6)	Value is outside the maximum and minimum limits
AT_ERR_NOTIMPLEMENTED (2)	Feature has not been implemented for the chosen camera
AT_ERR_READONLY (3)	Feature is read only
AT_ERR_NOTWRITABLE (5)	Feature is currently not writable
AT_ERR_NOTREADABLE (4)	Feature is currently not readable
AT_ERR_NULL_FEATURE (20)	NULL feature name passed to function
AT_ERR_NULL_READABLE_VAR (23)	Readable not set
AT_ERR_NULL_WRITABLE_VAR (25)	Writable not set
AT_ERR_NULL_ISAVAILABLE_VAR (31)	Available not set
AT_ERR_NULL_VALUE (28)	NULL value returned from function
AT_ERR_NULL_MINVALUE (26)	NULL min value
AT_ERR_NULL_MAXVALUE (27)	NULL max value
AT_ERR_INVALIDHANDLE (12)	Invalid device handle passed to function
AT_ERR_NOMEMORY (37)	No memory has been allocated for the current action
AT_ERR_COMM (17)	An error has occurred while communicating with hardware

Float Feature	Description
AT_SUCCESS (0)	Function call has been successful
AT_ERR_OUTOFRANGE (6)	Value is outside the maximum and minimum limits
AT_ERR_NOTIMPLEMENTED (2)	Feature has not been implemented for the chosen camera
AT_ERR_READONLY (3)	Feature is read only
AT_ERR_NOTWRITABLE (5)	Feature is currently not writable
AT_ERR_NOTREADABLE (4)	Feature is currently not readable
AT_ERR_NULL_FEATURE (20)	NULL feature name passed to function
AT_ERR_NULL_READABLE_VAR (23)	Readable not set
AT_ERR_NULL_WRITABLE_VAR (25)	Writable not set
AT_ERR_NULL_ISAVAILABLE_VAR (31)	Available not set
AT_ERR_NULL_VALUE (28)	NULL value returned from function
AT_ERR_NULL_MINVALUE (26)	NULL min value
AT_ERR_NULL_MAXVALUE (27)	NULL max value
AT_ERR_INVALIDHANDLE (12)	Invalid device handle passed to function
AT_ERR_NOMEMORY (37)	No memory has been allocated for the current action
AT_ERR_COMM (17)	An error has occurred while communicating with hardware

Boolean Feature	Description
AT_SUCCESS (0)	Function call has been successful
AT_ERR_OUTOFRANGE (6)	The value passed to the function was not a valid boolean value i.e. 0 or 1.
AT_ERR_NOTIMPLEMENTED (2)	Feature has not been implemented for the chosen camera
AT_ERR_READONLY (3)	Feature is read only
AT_ERR_NOTWRITABLE (5)	Feature is currently not writable
AT_ERR_NOTREADABLE (4)	Feature is currently not readable
AT_ERR_NULL_FEATURE (20)	NULL feature name passed to function
AT_ERR_NULL_READABLE_VAR (23)	Readable not set
AT_ERR_NULL_WRITABLE_VAR (25)	Writable not set
AT_ERR_NULL_ISAVAILABLE_VAR (31)	Available not set
AT_ERR_NULL_VALUE (28)	NULL value returned from function
AT_ERR_INVALIDHANDLE (12)	Invalid device handle passed to function
AT_ERR_NOMEMORY (37)	No memory has been allocated for the current action
AT_ERR_COMM (17)	An error has occurred while communicating with hardware

Enumerated Feature	Description
AT_SUCCESS (0)	Function call has been successful
AT_ERR_OUTOFRANGE (6)	The index passed to the function was either less than zero or greater than or equal to the number of implemented options.
AT_ERR_NOTIMPLEMENTED (2)	Feature has not been implemented for the chosen camera
AT_ERR_READONLY (3)	Feature is read only
AT_ERR_NOTWRITABLE (5)	Feature is currently not writable
AT_ERR_NOTREADABLE (4)	Feature is currently not readable
AT_ERR_INDEXNOTAVAILABLE (7)	Index is currently not available
AT_ERR_INDEXNOTIMPLEMENTED (8)	Index is not implemented for the chosen camera
AT_ERR_STRINGNOTAVAILABLE (18)	Index / String is not available
AT_ERR_STRINGNOTIMPLEMENTED (19)	Index / String is not implemented for the chosen camera
AT_ERR_NULL_FEATURE (20)	NULL feature name passed to function
AT_ERR_NULL_READABLE_VAR (23)	Readable not set
AT_ERR_NULL_WRITABLE_VAR (25)	Writable not set
AT_ERR_NULL_ISAVAILABLE_VAR (31)	Available not set
AT_ERR_NULL_VALUE (28)	NULL value returned from function
AT_ERR_NULL_COUNT_VAR (30)	NULL feature count
AT_ERR_NULL_IMPLEMENTED_VAR (22)	Feature not implemented
AT_ERR_INVALIDHANDLE (12)	Invalid device handle passed to function
AT_ERR_NOMEMORY (37)	No memory has been allocated for the current action
AT_ERR_COMM (17)	An error has occurred while communicating with hardware

Command Feature	Description
AT_SUCCESS (0)	Function call has been successful
AT_ERR_NOTIMPLEMENTED (2)	Feature has not been implemented for the chosen camera
AT_ERR_NOTWRITABLE (5)	Feature is currently not executable
AT_ERR_NULL_FEATURE (20)	NULL feature name passed to function
AT_ERR_NULL_READABLE_VAR (23)	Readable not set
AT_ERR_NULL_WRITABLE_VAR (25)	Writable not set
AT_ERR_NULL_ISAVAILABLE_VAR (31)	Available not set
AT_ERR_NULL_VALUE (28)	NULL value returned from function
AT_ERR_INVALIDHANDLE (12)	Invalid device handle passed to function
AT_ERR_NOMEMORY (37)	No memory has been allocated for the current action
AT_ERR_COMM (17)	An error has occurred while communicating with hardware

Buffer Control	Description
AT_SUCCESS (0)	Function call has been successful
AT_ERR_TIMEDOUT (13)	The AT_WaitBuffer function timed out while waiting for data arrive in output queue
AT_ERR_BUFFERFULL (14)	The input queue has reached its capacity
AT_ERR_INVALIDSIZE (15)	The size of a queued buffer did not match the frame size
AT_ERR_INVALIDALIGNMENT (16)	A queued buffer was not aligned on an 8-byte boundary
AT_ERR_HARDWARE_OVERFLOW (100)	The software was not able to retrieve data from the card or camera fast enough to avoid the internal hardware buffer bursting.
AT_ERR_NOMEMORY (37)	No memory has been allocated for the current action
AT_ERR_NODATA (11)	No Internal Event or Internal Error
AT_ERR_COMM (17)	An error has occurred while communicating with hardware
AT_ERR_NULL_QUEUE_PTR (34)	Pointer to queue is NULL
AT_ERR_NULL_WAIT_PTR (35)	Wait pointer is NULL
AT_ERR_NULL_PTRSIZE (36)	Pointer size is NULL

Feature Callback	Description
AT_SUCCESS (0)	Function call has been successful
AT_ERR_NULL_FEATURE (20)	NULL feature name passed to function
AT_ERR_NULL_EVCALLBACK (33)	EvCallBack parameter is NULL
AT_ERR_NOTIMPLEMENTED (2)	Feature has not been implemented for the chosen camera
AT_ERR_INVALIDHANDLE (12)	The size of a queued buffer did not match the frame size

SECTION 4 FEATURES

4.1 CAMERA SUPPORT

SDK3 currently supports only the Andor Neo and Zyla cameras, and the Apogee family of cameras. The features that are available for these cameras are outlined below. There is also a software module called **SimCam** that simulates the functionality of a camera. The SimCam module can be useful to prototype an application where availability of a real camera may be limited. To use SimCam, you should copy the atdevsimcam.dll file into your application directory. On initialisation of the SDK there will be two SimCam cameras available. The features available under SimCam are also outlined below.

4.2 FEATURE REFERENCE

Feature	Type	Description	Availability
AccumulateCount	Integer	Sets the number of images that should be summed to obtain each image in sequence.	Neo Zyla iStar-sCMOS
AcquiredCount	Integer	Dynamically incrementing count during an image sequence.	Apogee
AcquisitionStart	Command	Starts an acquisition.	SimCam, Neo Zyla, iStar-sCMOS
AcquisitionStop	Command	Stops an acquisition.	SimCam, Neo Zyla, iStar-sCMOS
AlternatingReadoutDirection	Boolean	Configures whether the sensor will be read out in alternating directions. See SensorReadDirection .	Zyla
AOIBinning	Enumerated	Sets up pixel binning on the camera. <u>Options:</u> <ul style="list-style-type: none"> • 1x1 • 2x2 • 3x3 • 4x4 • 8x8 See Section 4.6 Area of Interest .	Neo Zyla iStar-sCMOS
AOIHBin	Integer	Configures the Horizontal Binning of the sensor area of interest. See Section 4.6 Area of Interest .	SimCam Apogee, Zyla, iStar-sCMOS
AOIHeight	Integer	Configures the Height of the sensor area of interest in super-pixels. See Section 4.6 Area of Interest .	SimCam, Neo Zyla, Apogee, iStar-sCMOS
AOILayout	Enumerated	<u>Options:</u> <ul style="list-style-type: none"> • Image • Kinetics • TDI • Multitrack 	Apogee, Zyla, iStar-sCMOS
AOILeft	Integer	Configures the left hand coordinate of the sensor area of interest in sensor pixels. See Section 4.6 Area of Interest .	SimCam, Neo Zyla, Apogee, iStar-sCMOS
AOIStride	Integer	The size of one row in the image in bytes. Extra padding bytes may be added to the end of each line after pixel data to comply with line size granularity restrictions imposed by the underlying hardware interface. See Section 4.3 Image Format .	Neo Zyla, iStar-sCMOS
AOITop	Integer	Configures the top coordinate of the sensor area of interest in sensor pixels. See Section 4.6 Area of Interest .	SimCam, Neo Zyla, Apogee, iStar-sCMOS
AOIVBin	Integer	Configures the Vertical Binning of the sensor area of interest. See Section 4.6 Area of Interest .	SimCam Apogee, Zyla, iStar-sCMOS
AOIWidth	Integer	Configures the Width of the sensor area of interest in super-pixels.	SimCam, Neo Zyla, Apogee,

Feature	Type	Description	Availability
		See Section 4.6 Area of Interest .	iStar-sCMOS
AuxiliaryOutSource	Enumerated	Configures which signal appears on the auxiliary output pin. <u>Options:</u> <ul style="list-style-type: none">• FireRow1• FireRowN• FireAll• FireAny	Neo Zyla
AuxOutSourceTwo	Enumerated	AuxOutSourceTwo is a configurable output available to the user on the D-type (Zyla) and SMA (iStar-sCMOS). <u>Options:</u> <ul style="list-style-type: none">• ExternalShutterControl• FrameClock• RowClock• ExposedRowClock	Zyla, iStar-sCMOS
BackoffTemperatureOffset	Float	The Backoff temperature offset of the cooler subsystem.	Apogee
Baseline	Integer	Returns the baseline level of the image with current settings	Neo Zyla iStar-sCMOS
BitDepth	Enumerated	Returns the number bits used to store information about each pixel of the image. . Supported Bit Depth will be dependent on the camera. <u>Options Neo/Zyla:</u> <ul style="list-style-type: none">• 11 Bit or 12 Bit• 16 Bit <u>Options Apogee:</u> <ul style="list-style-type: none">• 12 Bit (Not available AltaF/Aspen/Ascent)• 16 Bit For AltaU/E this is determined by PixelReadoutRate: <ul style="list-style-type: none">• "Normal" -> "16-bit"• "Fast" -> "12-bit"	Neo Zyla, Apogee, iStar-sCMOS
BufferOverflowEvent	Integer	When enabled this will notify the user that the image buffer on the camera has been exceeded, causing the current acquisition to stop.	Neo Zyla iStar-sCMOS
BytesPerPixel	Floating Point	Returns the calculated bytes per pixel. This is read only.	Neo Zyla iStar-sCMOS
CameraAcquiring	Boolean	Returns whether or not an acquisition is currently acquiring.	SimCam, Neo Zyla, iStar-sCMOS
CameraDump	Command	Dumps current hardware configuration information to file in the executable directory. File is called camdump-[Serial Number]	Neo Zyla iStar-sCMOS
CameraFamily	String	Returns the family of the camera.	Apogee
CameraMemory	Integer	Returns the amount of available memory for storing images, in bytes.	Apogee
CameraModel	String	Returns the camera model.	SimCam, Neo Zyla, iStar-sCMOS
CameraName	String	Returns the name of the camera.	Neo, Zyla, Apogee, iStar-sCMOS
CameraPresent	Boolean	Returns whether the camera is connected to the system. Register a callback to this feature to be notified if the camera is disconnected.	
ColourFilter	Enumerated	Controls which colour filter is enabled. <u>Options:</u> <ul style="list-style-type: none">• None	Apogee

Feature	Type	Description	Availability
		<ul style="list-style-type: none"> • Bayer • TRUESENSE 	
ControllerID	String	Returns a unique identifier for the camera controller device. i.e. Frame grabber over Cameralink	Neo (CL) Zyla (CL)
CoolerPower	Double	Percentage of maximum power being used by the cooler.	Apogee
CycleMode	Enumerated	<p>Configures whether the camera will acquire a fixed length sequence or a continuous sequence. In Fixed mode the camera will acquire 'FrameCount' number of images and then stop automatically. In Continuous mode the camera will continue to acquire images indefinitely until the 'AcquisitionStop' command is issued.</p> <p>Options:</p> <ul style="list-style-type: none"> • Fixed • Continuous 	SimCam, Neo Zyla, iStar-sCMOS
DDGIOCEnable	Boolean	Enable or disable the integrate on chip option.	iStar-sCMOS
DDGIOCNumberOfPulses	Integer	Configures the number of pulses generated when the integrate on chip option is enabled.	iStar-sCMOS
DDGIOCPeiod	Integer	Configures the period (in picoseconds) of the pulses generated when the integrate on chip option is enabled.	iStar-sCMOS
DDGOutputDelay	Integer	Configures the delay (in picoseconds) for the currently selected DDG output.	iStar-sCMOS
DDGOutputEnable	Boolean	Enable or disable the currently selected DDG output.	iStar-sCMOS
DDGOutputStepEnable	Boolean	Enable or disable DDG step mode for the currently selected DDG output.	iStar-sCMOS
DDGOpticalWidthEnable	Boolean	Enable or disable the optical gate width option.	iStar-sCMOS
DDGOutputPolarity	Enumerated	<p>Configures the polarity for the currently selected DDG output.</p> <p>Options:</p> <ul style="list-style-type: none"> • Negative • Positive 	iStar-sCMOS
DDGOutputSelector	Enumerated	<p>Selects the DDG output that you subsequently wish to configure using the DDG output Features.</p> <p>Options:</p> <ul style="list-style-type: none"> • Gater • Output A • Output B • Output C 	iStar-sCMOS
DDGOutputWidth	Integer	Configures the width (in picoseconds) for the currently selected DDG output.	iStar-sCMOS
DDGStepCount	Integer	Configures the number of frames that will use the same step value. If you are using accumulations and you require each of the accumulated frames to use the same step values then this feature should be set to the number of accumulations.	iStar-sCMOS
DDGStepDelayCoefficientA	Floating Point	Configures coefficient A for the delay DDG step.	iStar-sCMOS
DDGStepDelayCoefficientB	Floating Point	Configures coefficient B for the delay DDG step.	iStar-sCMOS
DDGStepDelayMode	Enumerated	<p>Configures the delay DDG step mode.</p> <p>Options:</p> <ul style="list-style-type: none"> • Off • Constant • Exponential • Logarithmic 	iStar-sCMOS
DDGStepEnabled	Boolean	Indicates when DDG step is enabled.	iStar-sCMOS
DDGStepUploadProgress	Integer	Indicates the current progress of the DDG step upload.	iStar-sCMOS
DDGStepUploadRequired	Boolean	Indicates if a DDG step upload is required.	iStar-sCMOS
DDGStepWidthCoefficientA	Floating	Configures coefficient A for the width DDG step.	iStar-sCMOS

Feature	Type	Description	Availability
	Point		
DDGStepWidthCoefficientB	Floating Point	Configures coefficient B for the width DDG step.	iStar-sCMOS
DDGStepWidthMode	Enumerated	Configures the width DDG step mode. <u>Options:</u> <ul style="list-style-type: none"> • Off • Constant • Exponential • Logarithmic 	iStar-sCMOS
DDGStepUploadModeValues	Command	Uploads the DDG step values to the DDG frame buffer.	iStar-sCMOS
DDR2Type	String	Returns the amount of DDR2 Memory installed.	Apogee
DeviceCount	Integer	Returns the number of cameras detected.	System
DeviceVideoIndex	Integer	Returns the /dev/videoN number in Linux only	Neo (CL) Zyla (CL)
DisableShutter	Boolean	Enable or Disable shutter during acquisition. Overridden by ForceShutterOpen.	Apogee
DriverVersion	String	Returns USB driver version if USB Interface. Returns cURL library version if Ethernet interface.	Apogee
ElectronicShutteringMode	Enumerated	Configures which on-sensor electronic shuttering mode is used. <ul style="list-style-type: none"> • For pulsed or fast moving images <i>Global</i> shuttering is recommended. • For the highest frame rates and best noise performance <i>Rolling</i> is recommended. <u>Options:</u> <ul style="list-style-type: none"> • Rolling (not available in iStar-sCMOS) • Global 	SimCam, Neo Zyla, iStar-sCMOS
EventEnable	Boolean	Enable or Disable the acquisition event selected via the EventSelector feature.	Neo Zyla iStar-sCMOS
EventsMissedEvent	Integer	When enabled this will notify the user that an acquisition event, which the user registered a callback for, has been missed.	Neo Zyla iStar-sCMOS
EventSelector	Enumerated	Selects the acquisition events you wish to enable or disable using the EventEnable feature. <u>Options:</u> <ul style="list-style-type: none"> • ExposureEndEvent • ExposureStartEvent • RowNExposureEndEvent (not available in iStar-sCMOS) • RowNExposureStartEvent (not available in iStar-sCMOS) • EventsMissedEvent • BufferOverflowEvent 	Neo Zyla iStar-sCMOS
ExposedPixelHeight	Integer	Configures the exposure window in pixels.	Zyla iStar-sCMOS
ExposureTime	Floating Point	The requested exposure time in seconds. <i>Note: In some modes the exposure time can also be modified while the acquisition is running.</i>	SimCam, Neo Zyla, iStar-sCMOS
ExposureEndEvent	Integer	When enabled this will notify the user on the Negative edge of the FIRE in Global Shutter and FIRE of Row 1 in Rolling Shutter.	Neo Zyla iStar-sCMOS
ExposureStartEvent	Integer	When enabled this will notify the user on the Positive edge of the FIRE in Global Shutter and FIRE of Row 1 in Rolling Shutter.	Neo Zyla iStar-sCMOS
ExternalIOReadout	Boolean	When TRUE, the readout of the camera is no longer started by the external shutter. Instead, Pin 5 "External Readout Start" is used to start the readout. The default value of this variable after initialization is FALSE.	Apogee

Feature	Type	Description	Availability
		Not available on Ascent. Must use IOControl to cause the camera to consider Pin 5 to be "External Readout Start" and not user-defined.	
ExternalTriggerDelay	Floating Point	Sets the delay time between the camera receiving an external trigger and the acquisition start.	Zyla iStar-sCMOS
FanSpeed	Enumerated	Configures the speed of the fan in the camera. <u>Options SimCam/Neo/Zyla/iStar-sCMOS:</u> <ul style="list-style-type: none"> • Off • Low (Neo & SimCam Only) • On <u>Options Apogee:</u> <ul style="list-style-type: none"> • Off • Low • Medium • High 	SimCam, Neo Zyla, Apogee, iStar-sCMOS
FastAOIFrameRateEnable	Boolean	Enables faster framerates at small AOIs.	Zyla, Neo, iStar-sCMOS
FirmwareVersion	String	Returns the camera firmware version	Neo Zyla Apogee iStar-sCMOS
ForceShutterOpen	Boolean	Choose whether to force the shutter to open. Overrides DisableShutter.	Apogee
FrameCount	Integer	Configures the number of images to acquire in the sequence. The value of FrameCount must be any value which is a multiple of AccumulateCount. This ensures the accumulation contains the correct number of frames. When this feature is unavailable then the camera does not currently support fixed length series, therefore you must explicitly abort the acquisition once you have acquired the amount of frames required.	SimCam Neo Zyla, Apogee, iStar-sCMOS
FrameInterval	Floating Point	The interval in seconds between the end of readout of one image to the beginning of exposure of the next.	Apogee
FrameIntervalTiming	Boolean	Configures whether the timing of image acquisition is determined by FrameInterval or FrameRate	Apogee
FrameRate	Floating Point	Configures the frame rate in Hz at which each image is acquired during any acquisition sequence. This is the rate at which frames are acquired by the camera which may be different from the rate at which frames are delivered to the user. For example when AccumulateCount has a value other than 1, the apparent frame rate will decrease proportionally.	SimCam, Neo Zyla, Apogee, iStar-sCMOS
FullAOIControl	Boolean	Indicates whether or not the camera supports arbitrary AOI selection. If this feature is false then the supported AOI's are limited to those listed in Section 4.6 Area of Interest .	Neo Zyla iStar-sCMOS
GateMode	Enumerated	Selects the photocathode gating mode. <u>Options:</u> <ul style="list-style-type: none"> • CWOn • CWOff • FireOnly • GateOnly • FireAndGate • DDG 	iStar-sCMOS
HeatSinkTemperature	Float	Returns the current Heatsink Temperature. Not available for Ascent.	Apogee
I2CAddress	Integer	Configures the address of the I ² C device.	iStar-sCMOS
I2CByte	Integer	Get/Set the I ² C byte currently selected by the	iStar-sCMOS

Feature	Type	Description	Availability
		I2CByteSelector.	
I2CByteCount	Integer	Configures the number of I ² C bytes to read/write.	iStar-sCMOS
I2CByteSelector	Integer	Selects the I ² C byte that you subsequently wish to Get/Set using the I2CByte Feature.	iStar-sCMOS
I2CRead	Command	Performs the configured I ² C read.	iStar-sCMOS
I2CWrite	Command	Performs the configured I ² C write.	iStar-sCMOS
ImageSizeBytes	Integer	Returns the buffer size in bytes required to store the data for one frame. This will be affected by the Area of Interest size, binning and whether metadata is appended to the data stream.	SimCam, Neo Zyla, iStar-sCMOS
InputVoltage	Float	Returns the operating input voltage to the camera.	Apogee
InsertionDelay	Enumerated	Configures the length of the insertion delay. <u>Options:</u> <ul style="list-style-type: none"> • Normal • Fast 	iStar-sCMOS
InterfaceType	String	Returns the camera interface type. Current types: <u>Options Neo/Zyla:</u> <ul style="list-style-type: none"> • USB3 • CL 3 Tap • CL 2x5 Tap • CL 10 Tap <u>Options iStar-sCMOS:</u> <ul style="list-style-type: none"> • USB3 <u>Options Apogee:</u> <ul style="list-style-type: none"> • USB2 • Ethernet 	Neo Zyla Apogee iStar-sCMOS
IOControl	Enum	Configures whether selected IO is default or user defined. <u>Options</u> <ul style="list-style-type: none"> • Default • User 	Apogee
IODirection	Boolean or Enumerated	Configures whether selected IO is input or output. Cannot be changed for AltaF or Ascent. <u>Options</u> <ul style="list-style-type: none"> • Input / 0 • Output / 1 	Apogee
IOState	Boolean	Configures whether selected IO is enabled or disabled.	Apogee
IOInvert	Boolean	Indicates whether or not the operation of the IO Pin selected through the <i>IO Selector</i> Feature is inverted.	Neo Zyla iStar-sCMOS
IOSelector	Enumerated	Selects the IO Pin that you subsequently wish to configure using the <i>IO Invert</i> Feature. <u>Options Neo/Zyla:</u> <ul style="list-style-type: none"> • Fire 1 • Fire N (Zyla Only) • Aux Out 1 • Arm • External Trigger • Fire N and 1 (deprecated) <u>Options iStar-sCMOS:</u> <ul style="list-style-type: none"> • Aux Out 2 • External Trigger 	Neo Zyla Apogee iStar-sCMOS

Feature	Type	Description	Availability
		Options Apogee: <ul style="list-style-type: none"> External Trigger (I/O Signal 1: TriggerNormal) Fire (I/O Signal 2: Shutter Output) Shutter Strobe (I/O Signal 3: Shutter Strobe) External Exposure (I/O Signal 4 ExternalShutter) External Readout (I/O Signal 5: ExternalIOReadout) 	
IRPreFlashEnable	Boolean	When TRUE, the camera normalizes the sensor before an image is taken with a flash of IR. Not available for Ascent or other Interline transfer CCDs.	Apogee
KeepCleanEnable	Boolean	Enables/Disables any flushing command sent by the driver.	Apogee
KeepCleanPostExposure Enable	Boolean	Enables/Disables the camera control firmware to/from immediately beginning an internal flush cycle after an exposure.	Apogee
LineScanSpeed	Floating Point	Configures the number of rows read per second.	Zyla
LUTIndex	Integer	Sets the position in the LUT to read/write a new pixel map	
LUTValue	Integer	Sets the value in LUT in position specified by "LUT Index"	
MaxInterfaceTransferRate	Float	Returns the maximum sustainable transfer rate of the interface for the current shutter mode and AOI.	Neo Zyla iStar-sCMOS
MCPGain	Integer	Controls the voltage applied across the microchannel plate.	iStar-sCMOS
MCPIntelligate	Boolean	Enables or disables the MCP Intelligate mode. Only available when normal insertion delay is selected.	iStar-sCMOS
MCPVoltage	Integer	Indicates the current MCP voltage.	iStar-sCMOS
MetadataEnable	Boolean	Enable metadata. This is a global flag which will enable inclusion of metadata in the data stream as described in Section 4.5 Metadata . When this flag is enabled the data stream will always contain the MetadataFrame information. This will override the subsequent metadata settings when disabled. For example: If this feature is disabled and MetadataTimestamp is enabled, then metadata will not be included in the data stream. For example: If this feature is enabled and MetadataTimestamp is disabled, then metadata will be included in the data stream, but without timestamp information.	Neo Zyla iStar-sCMOS
MetadataFrame	Boolean	Indicates whether the MetadataFrame information is included in the data stream. This is read only and is automatically sent if metadata is enabled.	Neo Zyla iStar-sCMOS
MetadataTimestamp	Boolean	Enables inclusion of timestamp information in the metadata stream. The timestamp indicates the time at which the exposure for the frame started.	Neo Zyla iStar-sCMOS
MicrocodeVersion	String	Returns a revision code for the internal USB firmware within the camera head. Only available on USB interface.	Apogee
MultitrackBinned	Boolean	Configures whether the currently selected multitrack will be binned or not. Default state is set to true.	Zyla iStar-sCMOS
MultitrackCount	Integer	It is in the range 1-256. When set to 0, multitrack is disabled.	Zyla iStar-sCMOS
MultitrackEnd	Integer	Configures the row at which the currently selected multitrack ends.	Zyla iStar-sCMOS
MultitrackSelector	Integer	Selects multitrack index. It is in the range 0-255.	Zyla iStar-sCMOS
MultitrackStart	Integer	Configures the row at which the currently selected multitrack begins.	Zyla iStar-sCMOS
Overlap	Boolean	Enables overlap readout mode.	Neo Zyla Apogee

Feature	Type	Description	Availability
			iStar-sCMOS
PIVEnable	Boolean	Enables PIV mode.	iStar-sCMOS
PixelCorrection	Enumerated	Configures the pixel correction to be applied. <u>Options:</u> • Raw	SimCam
PixelEncoding	Enumerated	Configures the format of data stream. See Section 4.4 Pixel Encoding . <u>Neo, Zyla, iStar-sCMOS and SimCam Options:</u> Mono12 Mono12Packed Mono16 Mono32 <u>SimCam only Options:</u> Mono8 RGB8Packed Mono12Coded Mono12CodedPacked Mono22Parallel Mono22PackedParallel See Section 4.7 PixelEncoding and PreAmpGainControl for the dependency between this feature and the PreAmpGainControl feature.	SimCam, Neo, Zyla, iStar-sCMOS
PixelHeight	Floating Point	Returns the height of each pixel in micrometers.	SimCam, Neo, Zyla, Apogee, iStar-sCMOS
PixelReadoutRate	Enumerated	Configures the rate of pixel readout from the sensor. <u>Options sCMOS:</u> • 280 MHz • 200 MHz (Neo Only - deprecated) • 100 MHz <u>Options SimCam:</u> • 550 MHz <u>Options Apogee:</u> • Normal • Fast (Not Available on AltaE)	SimCam, Neo, Zyla, Apogee, iStar-sCMOS
PixelWidth	Floating Point	Returns the width of each pixel in micrometers.	SimCam, Neo, Zyla, Apogee, iStar-sCMOS
PortSelector	Integer	Configures the port being used. Not implemented if fixed Single readout. Min is always 0. Max is {0, 1, 3} depending on {Single, Dual, Quad} readout. If the camera is switchable between Single and Dual readout then the value of PortSelector required to set the gain/offset for Single readout is camera-dependent (i.e. it could be 0 or 1).	Apogee
PreAmpGain	Enumerated	Configures the gain applied to the gain channel selected through the <i>Pre Amp Gain Selector</i> Feature. <u>Options:</u> • x1 • x2 • x10 • x30	SimCam, Neo (deprecated)

Feature	Type	Description	Availability
PreAmpGainChannel	Enumerated	Configures which pre amp gain channel(s) will be used for reading out the sensor. <u>Options:</u> <ul style="list-style-type: none">• High• Low• Both	SimCam, Neo (deprecated)
PreAmpGainControl	Enumerated	Wrapper Feature to simplify access to the PreAmpGain, PreAmpGainChannel and PreAmpGainSelector features. See Section 4.7 PixelEncoding and PreAmpGainControl for the dependency between this feature and the PixelEncoding feature. This feature is deprecated and should be replaced by the SimplePreAmpGainControl feature as some of the options may not be supported. <u>Options:</u> <ul style="list-style-type: none">• Gain 1 (11 bit)• Gain 2 (11 bit)• Gain 3 (11 bit)• Gain 4 (11 bit)• Gain 1 Gain 3 (16 bit)• Gain 1 Gain 4 (16 bit)• Gain 2 Gain 3 (16 bit)• Gain 2 Gain 4 (16 bit)	Neo (deprecated)
PreAmpGainValue	Integer	This is the value of the pre-amplifier gain, for the currently-selected ADC/channel pair as selected by PortSelector if PortSelector is implemented. Only effective for Fast PixelReadoutRate for AltaU/AltaE/AltaF/Aspen.	Apogee
PreAmpGainSelector	Enumerated	Selects the gain channel that you subsequently wish to configure using the <i>Pre Amp Gain</i> Feature. <u>Options:</u> <ul style="list-style-type: none">• High• Low	SimCam, Neo (deprecated)
PreAmpOffsetValue	Integer	This is the value of the pre-amplifier offset, for the currently-selected ADC/channel pair as selected by PortSelector, if PortSelector is implemented. Only effective for Fast PixelReadoutRate for AltaU/AltaE/AltaF/Aspen.	Apogee
PreTriggerEnable	Boolean	Enable or disable the pre-trigger input.	iStar-sCMOS
ReadoutTime	Floating Point	This feature will return the time to readout data from a sensor.	Neo Zyla iStar-sCMOS
RollingShutterGlobalClear	Boolean	Enables Rolling Shutter Global Clear readout mode.	Zyla
RowNExposureEndEvent	Integer	When enabled this will notify the user on the Negative edge of the FIRE of ROW N in Rolling Shutter.	Neo Zyla
RowNExposureStartEvent	Integer	When enabled this will notify the user on the Positive edge of the FIRE of ROW N in Rolling Shutter.	Neo Zyla
RowReadTime	Floating Point	Configures the time in seconds to read a single row.	Zyla
ScanSpeedControlEnable	Boolean	Configures whether the LineReadSpeed and RowReadTime can be altered.	Zyla
SensorCooling	Boolean	Configures the state of the sensor cooling. Cooling is disabled by default at power up and must be enabled for the camera to achieve its target temperature. The actual target temperature can be set with the TemperatureControl feature where available for example on the Neo camera.	SimCam, Neo Zyla, Apogee, iStar-sCMOS
SensorHeight	Integer	Returns the height of the sensor in pixels.	SimCam, Neo Zyla, Apogee,

Feature	Type	Description	Availability
			iStar-sCMOS
SensorModel	String	Returns the sensor model installed in the camera.	Apogee
SensorReadoutMode	Enumerated	Configures the direction in which the sensor will be read out. <u>Options</u> <ul style="list-style-type: none"> • Bottom Up Sequential • Bottom Up Simultaneous • Centre Out Simultaneous • Outside In Simultaneous • Top Down Sequential • Top Down Simultaneous 	Zyla
SensorType	Enumerated	Returns true for CCD, False for CMOS. <u>Options:</u> <ul style="list-style-type: none"> • CCD • CMOS 	Apogee
SensorTemperature	Floating Point	Read the current temperature of the sensor.	SimCam, Neo Zyla, Apogee, iStar-sCMOS
SensorWidth	Integer	Returns the width of the sensor in pixels.	SimCam, Neo Zyla, Apogee, iStar-sCMOS
SerialNumber	String	Returns the camera serial number.	SimCam, Neo Zyla, iStar-sCMOS
ShutterAmpControl	Boolean	Disables the CCD voltage while the shutter strobe is high.	Apogee
ShutterMode	Enumerated	Controls the behavior of the shutter. <u>Options:</u> <ul style="list-style-type: none"> • Open • Closed • Auto 	Apogee
ShutterOutputMode	Enumerated	Controls the mode the external trigger will run in. External Shutter signal can either be set to high (open) or low (closed). ShutterOutput can be triggered by setting AuxOutSourceTwo to ExternalShutterControl. <u>Options:</u> <ul style="list-style-type: none"> • Open • Closed 	Zyla iStar-sCMOS
ShutterState	Boolean	Returns whether shutter is opened or closed.	Apogee
ShutterStrobePeriod	Double	Sets the period of the shutter strobe on pin3. Must use IOControl to cause the camera to consider Pin 3 to be "Shutter Strobe Output" and not user-defined.	Apogee
ShutterStrobePosition	Double	Sets the delay from the time the exposure begins to the time the rising edge of the shutter strobe period appears on pin 3. Must use IOControl to cause the camera to consider Pin 3 to be "Shutter Strobe Output" and not user-defined.	Apogee
ShutterTransferTime	Floating Point	Sets the time at which the shutter will be opened before an exposure starts. Can be set between 0 and 100ms at 0.5ms steps. Only available whenever the ExternalTriggerMode is set to 'Automatic'	Zyla iStar-sCMOS
SimplePreAmpGainControl	Enumerated	Wrapper Feature to simplify selection of the sensitivity and dynamic range options. This feature should be used as a replacement for the PreAmpGainControl feature as some of the options in the PreAmpGainControl feature are not supported on all cameras. Supported Bit Depth will be dependent on the camera. See Section 4.7	Neo Zyla iStar-sCMOS

Feature	Type	Description	Availability
		<p>PixelEncoding and PreAmpGainControl for the dependency between this feature and the PixelEncoding feature.</p> <p><u>Options:</u></p> <ul style="list-style-type: none"> 11-bit (high well capacity) Or 12-bit (high well capacity) 11-bit (low noise) Or 12-bit (low noise) 16-bit (low noise & high well capacity) 	
SoftwareTrigger	Command	Generates a software trigger in the camera. Used to generate each frame on the camera whenever the trigger mode is set to <i>Software</i> .	Neo Zyla iStar-sCMOS
SoftwareVersion	String	Returns the version of the SDK.	System
SpuriousNoiseFilter	Boolean	Enables or Disables the Spurious Noise Filter	Neo Zyla iStar-sCMOS
StaticBlemishCorrection	Boolean	Enables or Disables Static Blemish Correction	Neo Zyla iStar-sCMOS
Synchronous Triggering	Boolean	Configures whether external triggers are synchronous with the read out of a sensor row. Asynchronous triggering may result in data corruption in the row being digitised when the triggers occurs.	SimCam
TargetSensor Temperature	Floating Point	Configures the temperature to which the sensor will be cooled. To be used for cameras with no correction data (-50->25). Otherwise <i>TemperatureControl</i> should be used.	SimCam, Neo (deprecated)
Temperature Control	Enumerated	Allows the user to set the target temperature of the sensor based on a list of valid temperatures.	Neo
Temperature Status	Enumerated	<p>Reports the current state of cooling towards the Target Sensor Temperature. [Read Only]</p> <p><u>Options Neo/Zyla/Apogee/iStar-sCMOS:</u></p> <ul style="list-style-type: none"> Cooler Off Stabilised Cooling <p><u>Neo/Zyla/iStar-sCMOS Only:</u></p> <ul style="list-style-type: none"> Drift Not Stabilised Fault <p><u>Apogee Only:</u></p> <ul style="list-style-type: none"> Backoff 	Neo Zyla, Apogee iStar-sCMOS
TimestampClock	Integer	Reports the current value of the camera's internal timestamp clock. This same clock is used to timestamp images as they are acquired when the MetadataTimestamp feature is enabled. The clock is reset to zero when the camera is powered on and then runs continuously at the frequency indicated by the TimestampClockFrequency feature. The clock is 64-bits wide.	Neo Zyla iStar-sCMOS
TimestampClock Frequency	Integer	Reports the frequency of the camera's internal timestamp clock in Hz.	Neo Zyla iStar-sCMOS
TimestampClock Reset	Command	Resets the camera's internal timestamp clock to zero. As soon as the reset is complete the clock will begin incrementing from zero at the rate given by the TimestampClockFrequency feature.	Neo Zyla iStar-sCMOS
TransmitFrames	Boolean	If false all image data will be concatenated and obtained as a single download at the end of the sequence.	Apogee

Feature	Type	Description	Availability
		If true each image (or image row in the case of TDI) from a sequence will be available for download after it is read out (digitised).	
TriggerMode	Enumerated	<p>Allows the user to configure the camera trigger mode at a high level. If the trigger mode is set to <i>Advanced</i> then the <i>Trigger Selector</i> and <i>Trigger Source</i> feature must also be set.</p> <p><u>Neo, Zyla, iStar-sCMOS, SimCam and Apogee Options:</u></p> <ul style="list-style-type: none"> • Internal • Software (Not available on Apogee) • External • External Start • External Exposure <p><u>SimCam only Options:</u></p> <ul style="list-style-type: none"> • Advanced 	SimCam, Neo Zyla, Apogee, iStar-sCMOS
UsbProductId	Integer	<p>Returns the USB Product ID associated with the camera system.</p> <p>Only available for USB interface.</p>	Apogee
UsbDeviceId	Integer	<p>Returns the USB Device ID associated with the camera system.</p> <p>Only available on USB interface.</p>	Apogee
VerticallyCentreAOI	Boolean	Vertically centres the AOI in the frame. With this enabled, AOITop will be disabled.	Neo, Zyla, iStar-sCMOS

4.3 IMAGE FORMAT

Images are presented to the application in the general format shown in **Figure 1**. Pixels are returned row by row starting from the top row and with the leftmost pixel being sent first in each row. The number of pixels in each row can be obtained from the **AOIWidth** feature and the number of rows in the image can be obtained from the **AOIHeight** feature.

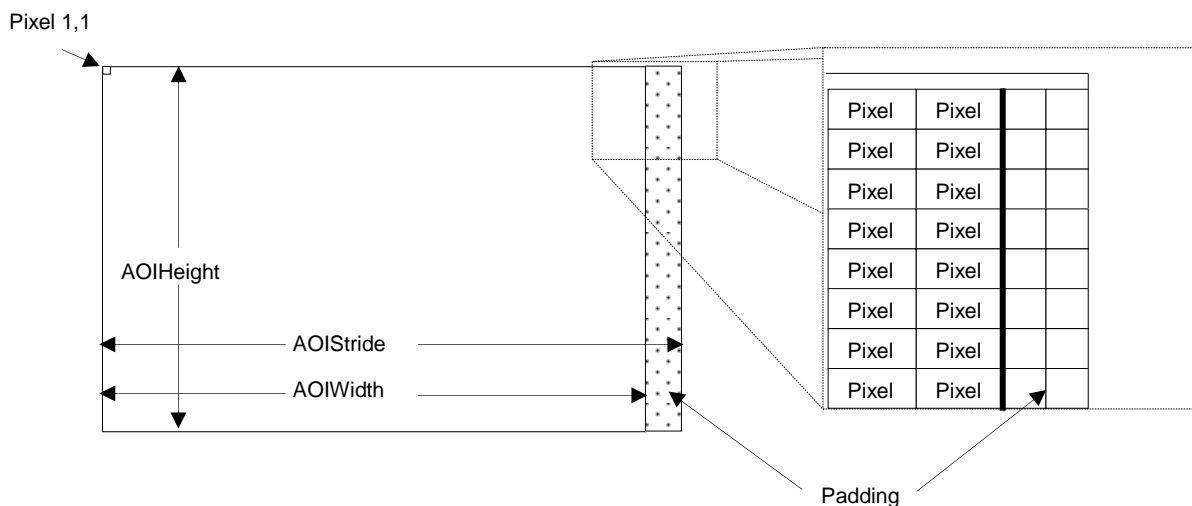


Figure 1: General Image Format

Stride

At the end of each row there may be additional padding bytes. This padding area does not contain any valid pixel data and should be skipped over when processing or displaying an image. This padding is necessary to ensure the image can be transferred over the interface between the camera and the PC and its size is dependent on the specific hardware interface being used as well as the current AOI settings. **Figure 2** shows what this padding looks like when viewing the raw data for an image in memory; in this example the pixels are 16-bit wide and the AOIWidth is 10.

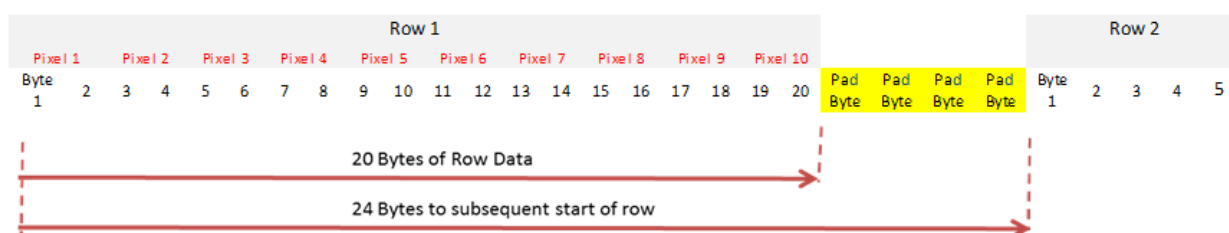


Figure 2: Padding between image rows in memory

To ensure that your application will operate successfully with any hardware interface or AOI configuration, you should make use of the **AOIStride** feature to skip over the padding at the end of each row. The AOIStride feature represents the total number of bytes that each row of the image contains and includes the memory necessary for pixel data plus any padding at the end of the row. The stride should be used during processing of an image to obtain the memory address of each row relative to the previous row. See the example code below. Note that AOIStride is measured in bytes whereas AOIWidth is measured in pixels.

```
//Get the next image from the SDK
AT_WaitBuffer(&ImageBuffer, &ImageSize, AT_INFINITE);

//Retrieve the dimensions of the image
AT_GetInt(Hndl, L"AOIStride", &Stride);
AT_GetInt(Hndl, L"AOIWidth", &Width);
```



```
AT_GetInt(Hndl, L"AOIHeight", &Height);

for (AT_64 Row=0; Row < Height; Row++) {

    //Cast the raw image buffer to a 16-bit array.
    //...Assumes the PixelEncoding is 16-bit.
    unsigned short* ImagePixels = reinterpret_cast<unsigned short*>(ImageBuffer);

    //Process each pixel in a row as normal
    for (AT_64 Pixel=0; Pixel < Width; Pixel++) {
        SomeProcessing(ImagePixels[Pixel]);
    }

    //Use Stride to get the memory location of the next row.
    ImageBuffer += Stride;
}
```

ATUtility Library

The **ATUtility library** provided with the SDK contains functionality that can be used to strip the padding from an image. Once this is done the image can be processed without concern for padding. Stripping the padding from an image will however incur some processing overhead. See example code below and **Section 6.1 ATUTILITY**.

```
//Get the next image from the SDK
AT_WaitBuffer(&ImageBuffer, &ImageSize, AT_INFINITE);

//Retrieve the dimensions of the image
AT_GetInt(Hndl, L"AOIStride", &Stride);
AT_GetInt(Hndl, L"AOIWidth", &Width);
AT_GetInt(Hndl, L"AOIHeight", &Height);

unsigned short ImagePixels[Width*Height];

//Use an atutility function to strip padding from the image
AT_ConvertBuffer(ImageBuffer, ImagePixels, Width, Height,
                Stride, L"Mono16", L"Mono16");

for (AT_64 Pixel=0; Pixel < Width*Height; Pixel++) {
    SomeProcessing(ImagePixels[Pixel]);
}
```

4.4 PIXEL ENCODING

There are several **Pixel Encoding** options available for the pixels in an image. Each of the formats is described below showing the pixel size, its layout in memory and sample C++ code for extracting pixel information out of the raw memory array.

In the descriptions below MSB refers to the most significant bits of the pixel, LSB refers to the least significant bits. ImageBuffer is the address of the start of the image in memory.

Mono12Packed

12-bit Monochrome Data, stored by packing two adjacent pixels into three bytes.

Bits 11:4 (MSB)		Bits 3:0 (LSB)	
ImageBuffer+0	Pixel A (MSB)		
ImageBuffer+1	Pixel B (LSB)	Pixel A (LSB)	
ImageBuffer+2	Pixel B (MSB)		
ImageBuffer+3	Pixel C (MSB)		
ImageBuffer+4	Pixel D (LSB)	Pixel C (LSB)	
ImageBuffer+5	Pixel D (MSB)		

```
PixelA = (ImageBuffer [0] << 4) + (ImageBuffer [1] & 0xF);
PixelB = (ImageBuffer [2] << 4) + (ImageBuffer [1] >> 4);
PixelC = (ImageBuffer [3] << 4) + (ImageBuffer [4] & 0xF);
PixelD = (ImageBuffer [5] << 4) + (ImageBuffer [4] >> 4);
```

Note: The atutility library can be used to easily convert an image in this format to Mono16. i.e.

```
unsigned short* ImagePixels [Width*Height];
AT_ConvertBuffer(ImageBuffer, ImagePixels, Width, Height,
                Stride, L"Mono12Packed", L"Mono16");
```

See **Section 6.1 ATUTILITY**.

Mono12

12-bit Monochrome Data, stored as 16-bit little-endian with zero padded upper bits.

Bits 11:8 (MSB)		Bits 7:0 (LSB)	
ImageBuffer+0	Pixel A (LSB)		
ImageBuffer+1	0	Pixel A (MSB)	
ImageBuffer+2	Pixel B (LSB)		
ImageBuffer+3	0	Pixel B (MSB)	
ImageBuffer+4	Pixel C (LSB)		
ImageBuffer+5	0	Pixel C (MSB)	

```
unsigned short* ImagePixels = reinterpret_cast<unsigned short*>(ImageBuffer);
PixelA = ImagePixels [0];
PixelB = ImagePixels [1];
PixelC = ImagePixels [2];
```

Mono16

16-bit Monochrome Data, stored as 16-bit little-endian.

<i>Bits 15:8 (MSB)</i>	<i>Bits 7:0 (LSB)</i>
ImageBuffer+0	Pixel A (LSB)
ImageBuffer+1	Pixel A (MSB)
ImageBuffer+2	Pixel B (LSB)
ImageBuffer+3	Pixel B (MSB)
ImageBuffer+4	Pixel C (LSB)
ImageBuffer+5	Pixel C (MSB)

```
unsigned short* ImagePixels = reinterpret_cast<unsigned short*>(ImageBuffer);
PixelA = ImagePixels [0];
PixelB = ImagePixels [1];
PixelC = ImagePixels [2];
```

Mono32

32-bit Monochrome Data, stored as 32-bit little-endian.

<i>Bits 31:24 (MSB)</i>	<i>Bits 23:16</i>	<i>Bits 15:8</i>	<i>Bits 7:0 (LSB)</i>
ImageBuffer+0	Pixel A (LSB)		
ImageBuffer+1	Pixel A		
ImageBuffer+2	Pixel A		
ImageBuffer+3	Pixel A (MSB)		
ImageBuffer+4	Pixel B (LSB)		
ImageBuffer+5	Pixel B		
ImageBuffer+6	Pixel B		
ImageBuffer+7	Pixel B (MSB)		

```
unsigned int* ImagePixels = reinterpret_cast<unsigned int*>(ImageBuffer);
PixelA = ImagePixels [0];
PixelB = ImagePixels [1];
```

Mono8 (Limited Availability)

8-bit Monochrome Data

<i>Bits 7:0</i>

ImageBuffer+0	Pixel A
ImageBuffer+1	Pixel B
ImageBuffer+2	Pixel C
ImageBuffer+3	Pixel D

```
unsigned char* ImagePixels = reinterpret_cast<unsigned char*>(ImageBuffer);
PixelA = ImagePixels [0];
PixelB = ImagePixels [1];
PixelC = ImagePixels [2];
ImagePixels [3];
```

4.5 METADATA

Metadata can be enabled through the **MetadataEnable** Boolean feature. When metadata is enabled extra information will be appended onto each image by the camera. The **ImageSizeBytes** feature will update to include the extra memory required for the metadata. By default metadata is disabled.

The features used to configure Metadata are:

- **MetadataEnable** – Enable inclusion of metadata information in the data stream.
- **MetadataFrameInfo** – Enable inclusion of frame information in the data stream.
- **MetadataTimestamp** – Enable inclusion of timestamp information in the data stream.
- **MetadataFrame** – Enable inclusion of image data in the data stream, this will always be included if metadata is enabled.

With Metadata enabled the format of the image stream is shown below.

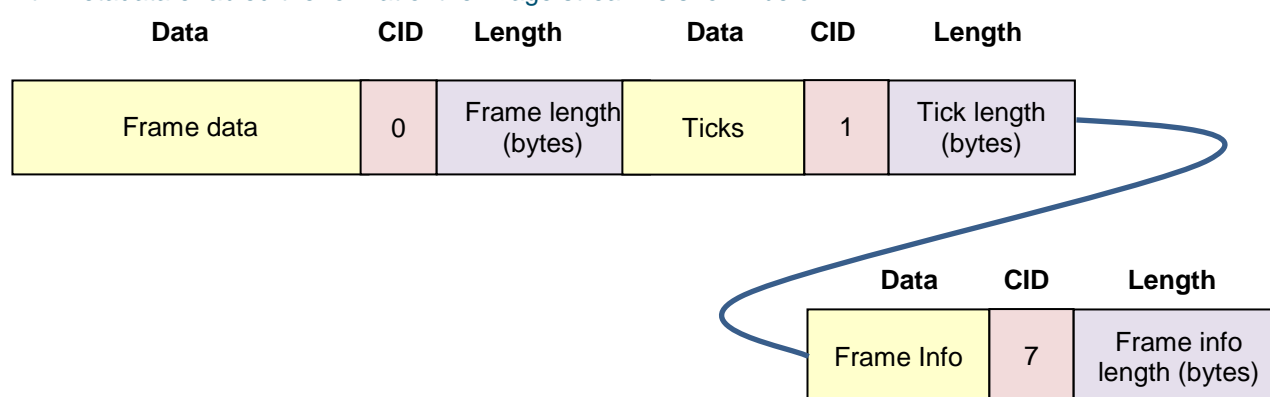


Figure 3: Format of the image stream with Metadata enabled.

Usage: Data, CID (Chunk ID), Length and then repeated as necessary.

Metadata is composed in blocks, with each block representing a particular type of metadata. Each block has an independent identifier called the **CID** (Chunk Identifier). As well as the CID, each block contains the actual metadata value and also a length field to facilitate parsing through the metadata blocks. The Length field in each block indicates the size of the metadata information plus the size of the CID. Note that the Length field is the last field in each metadata block. Parsing of metadata should be done in reverse, starting from the end of the data stream, and working back through the metadata until you reach the start of the data stream.

The three types of metadata block currently supported are **Frame (CID 0)**, **Timestamp (CID 1)** and **Frame Information (CID 7)**. The Frame metadata block simply contains the image data and is always enabled if metadata is enabled. The Timestamp metadata block contains a timestamp indicating the time at which the exposure for the frame was started. The Frame Information metadata block contains information on the structure of the image data.

All fields in the metadata blocks are stored little endian. i.e. least significant byte first.

The metadata format is described below:

Data

Data will vary depending on which CID it is. In this case it may be the actual image data or the timestamp information.

Chunk Identifier (CID)

A CID is used to label each block. Each CID is 4 bytes in length. The valid values for CID are shown below:

0 - Frame Data

- This represents the actual image.

1 – FPGA “Ticks”

- From camera power up, a 64 bit counter will count number of FPGA clocks or “Ticks”.
- The Ticks data will always be a 64 bit number – 8 bytes

7 – Frame Info

- Information about the frame including AOI, pixel encoding and stride.

Length

The length is a 4 byte number. This is where the length in bytes of the metadata block is stored. It includes the size of both the Data and the CID field.

Note

For the Frame Metadata block, Length is equal to the stride length x number of image rows, plus any padding at the end of the image plus the size of the CID field.

See the Tutorial section 2.1.9 Metadata for example code showing how to configure and parse metadata.

Timestamp Frequency

The frequency of the timestamp clock can be retrieved through the TimestampClockFrequency feature.

Timestamp Clock

The current value of the timestamp clock can be read directly from the camera by accessing the TimestampClock feature. This can be used to synchronise the timestamp attached to each image with an absolute calendar time. To do this the program should first read the current time from the PC clock, then read the TimestampClock feature. This reference point can then be used to find out the absolute time at which any image was acquired. Alternatively the TimestampClockReset feature can be executed which will reset the timestamp clock to zero.

Frame Info

The frame info block takes the form:

<i>Bits 63:48 (MSB)</i>	<i>Bits 47:32</i>	<i>Bits 31:24</i>	<i>Bits 23:16</i>	<i>Bits 15:0 (LSB)</i>
<i>AOI Height</i>	<i>AOI Width</i>	<i>0</i>	<i>Pixel Encoding</i>	<i>Stride</i>

AOI Height, AOI Width and Stride are all 16 bit numbers. Pixel encoding takes the values:

- 0 for Mono16
- 1 for Mono12
- 2 for Mono12Packed.
- 3 for Mono32

4.6 AREA OF INTEREST

The **Area of Interest (AOI)**, sometimes referred to as Region of Interest (ROI) is configured with the following features:

- AOIHBin, AOIVBin or AOIBinning
- AOIWidth
- AOILeft
- AOIHeight
- VerticallyCentreAOI
- AOITop

It is recommended that these features are configured in the order listed above as features towards the top of the list will override the values below them if the values are incompatible.

Super-Pixels

A **super-pixel** is the result of combining multiple sensor pixels into a single data pixel by binning the values from each sensor pixel together. The amount of binning in each direction is configured either by setting the AOIHBin and AOIVBin features or by using the AOIBinning feature. The AOIWidth and AOIHeight features are set and retrieved in units of super-pixels. Therefore, when binning is in use, the AOIWidth value will always indicate the number of data pixels that each row of the image data contains and not the number of pixels read off the sensor. The AOILeft and AOITop coordinates are specified in units of sensor pixels.

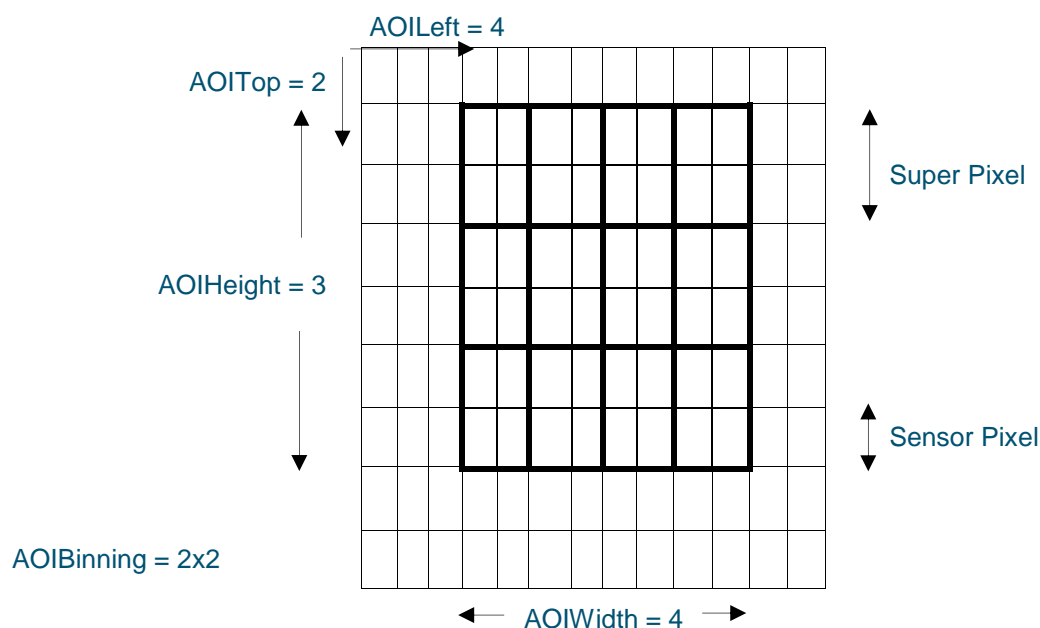


Figure 4: Configuring an AOI and presentation of Super-Pixels

Support for AOI Control

Some older versions of the Neo camera do not support full control over the AOI. In this case the FullAOIControl feature will return false.

If the FullAOIControl feature is not implemented or returns false for your camera then there are restrictions on the configurations of AOI's that can be used. The available AOI's are listed in the table below.

Width	Height	Top	Left(12bit)*	Left(16bit)**	X Centre (12 bit)	X Centre (16bit)
2592	2160	1	1	1	1297	1297
2544	2160	1	17	25	1289	1297
2064	2048	57	257	265	1289	1297
1776	1760	201	401	409	1289	1297
1920	1080	537	337	337	1297	1297
1392	1040	561	593	601	1289	1297
528	512	825	1025	1033	1289	1297
240	256	953	1169	1177	1289	1297
144	128	1017	1217	1225	1289	1297
2592	304	929	1	1	1297	1297

*12bit refers to mono12packed pixel encoding.

**16bit refers to mono12 or mono16 pixel encoding.

4.7 PIXELENCODING AND PREAMPAINCONTROL

When changing PreAmpGainControl, the PixelEncoding feature will automatically adjust to a default setting for that PreAmpGainControl. For example, if the PreAmpGainControl feature is changed from '12-bit (low noise)' to '16-bit (low noise & high well capacity)' then the PixelEncoding feature will switch automatically to Mono16. This will only occur if the PixelEncoding is currently set to Mono12 or Mono12Packed. Similarly, PixelEncoding will switch automatically to either Mono12 or Mono12Packed when the PreAmpGainControl feature is changed from '16-bit (low noise & high well capacity)' to '12-bit (low noise)'. When PixelEncoding is set to Mono32, it does not change when PreAmpGainControl is changed.

4.8 SENSOR COOLING

It is important to cool the temperature of the sCMOS sensor to reduce the amount of noise in the images captured, see example below of the same image at different sensor temperatures. Sensor cooling can be set with the Boolean feature SensorCooling. The sensor temperature can then be set with the Enumerated feature, TemperatureControl and read using the Float feature SensorTemperature. To check the status of the cooling mechanism, read the TemperatureStatus feature. The possible status options are:

- **Cooler Off** Cooling has been disabled.
- **Stabilised** Temperature has stabilised at Target Temperature.
- **Cooling** Temperature is approaching Target Temperature.
- **Drift** Temperature has drifted outside Target Range having stabilised.
- **Not Stabilised** Temperature is within Target Range but has not yet stabilised.
- **Fault** Temperature has been outside Target Range for a long period of time.

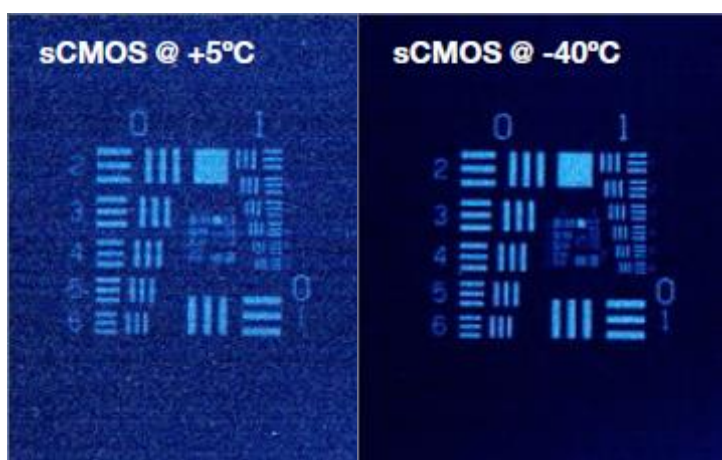


Figure 5: The effect of sensor cooling on noise: sCMOS sensor at +5°C and -40°C cooling.

```
double temperature = 0;
AT_SetBool(Hndl, L"SensorCooling", AT_TRUE);
AT_GetFloat(Hndl, L"SensorTemperature", &temperature);

cout << "Temperature: " << temperature << endl;

int temperatureCount = 0;
AT_GetEnumCount(Hndl, L"TemperatureControl", &temperatureCount);
AT_SetEnumIndex(Hndl, L"TemperatureControl", temperatureCount-1);

int temperatureStatusIndex = 0;
wchar_t* temperatureStatus[256];
AT_GetEnumIndex(Hndl, L"TemperatureStatus", &temperatureStatusIndex);
AT_GetEnumStringByIndex(Hndl, L"TemperatureStatus", temperatureStatusIndex,
temperatureStatus, 256);

while(wcscmp(L"Stabilised",temperatureStatus) != 0) {
    sleep(1);
    AT_GetEnumIndex(Hndl, L"TemperatureStatus", &temperatureStatusIndex);
    AT_GetEnumStringByIndex(Hndl, L"TemperatureStatus", temperatureStatusIndex,
temperatureStatus, 256);
    wcout << L"Temperature Status: " << temperatureStatus << endl;
}

cout << "Temperature Stabilised" << endl;
```


4.9 COMPARISON OF SDK2 AND SDK3

Action	SDK2	SDK3
Setting Exposure Time	<code>SetExposureTime (0.1);</code>	<code>AT_SetFloat (Hndl, L"ExposureTime", 0.1);</code>
Using External Trigger	<code>SetTriggerMode (1);</code>	<code>AT_SetEnumString (Hndl, L"TriggerMode", L"External");</code>
Starting an Acquisition	<code>StartAcquisition ();</code>	<code>AT_Command (Hndl, L"AcquisitionStart");</code>
Determining if Pre Amp Gain is supported	<pre>AndorCapabilities AndorCaps; GetCapabilties (&AndorCaps); if (AndorCaps.ulSetFunctions & AC_SETFUNCTION_PREAMPGAIN) Implemented = 1; else Implemented = 0;</pre>	<code>AT_IsImplemented (Hndl, L"PreAmpGain", &Implemented);</code>
Getting Current Exposure Time	<code>GetAcquisitionTimings (&Exposure, &AccCycleTime, &KinCycleTime);</code>	<code>AT_GetFloat (Hndl, L"ExposureTime", &Exposure);</code>
Setting Frame Rate (5fps)	<code>SetKineticCycleTime (0.2);</code>	<code>AT_SetFloat (Hndl, L"FrameRate", 5.0);</code>
Getting the Serial Number	<code>GetCameraSerialNumber (&Serial);</code>	<code>AT_GetString (Hndl, L"SerialNumber", Serial);</code>
Getting Current Trigger Mode	-- Not supported --	<code>AT_GetEnumIndex (Hndl, L"TriggerMode", Mode);</code>
AOI Setup	<code>SetImage (1, 1, 256, 384, 512, 768);</code>	<pre>AT_SetInt (Hndl, L"AOIHBin", 1); AT_SetInt (Hndl, L"AOIWidth", 128); AT_SetInt (Hndl, L"AOILeft", 256); AT_SetInt (Hndl, L"AOIVBin", 1); AT_SetInt (Hndl, L"AOIHeight", 256); AT_SetInt (Hndl, L"AOITop", 512);</pre>
Opening Second device	<code>GetCameraHandle (1, &SecondHandle)</code> <code>SetCurrentCamera (SecondHandle);</code>	<code>AT_Open (1, &SecondHandle);</code>
Getting Limits	<code>GetMinimumNumberInSeries (&Min)</code> -- Not supported --	<pre>AT_GetIntMin (Hndl, L"FrameCount", &Min); AT_GetIntMax (Hndl, L"FrameCount", &Max);</pre>

SECTION 5 FUNCTION REFERENCE

5.1 FUNCTION LISTING

This section provides a description of various reference functions available in SDK3.

5.1.1 AT_OPEN

```
int AT_Open(int DeviceIndex, AT_H* Handle)
```

Description

This function is used to open up a handle to a particular camera. The DeviceIndex parameter indicates the index of the camera that you wish to open and the handle to the camera is returned in the Handle parameter. This Handle parameter must be passed as the first parameter to all other functions to access the features or to acquire data from the camera.

5.1.2 AT_CLOSE

```
int AT_Close(AT_H Handle)
```

Description

This function is used to close a previously opened handle to a camera. The Handle parameter is the handle that was returned from the AT_Open function. The function should be called when you no longer wish to access the camera from your application usually at shutdown.

5.1.3 AT_ISIMPLEMENTED

```
int AT_IsImplemented(AT_H Hndl, AT_WC* Feature, AT_BOOL* Implemented)
```

Description

This function can be used to determine whether the camera has implemented the feature specified by the Feature parameter. On return the Implemented parameter will contain the value AT_FALSE or AT_TRUE. In the case that the feature is implemented the value of Implemented will be AT_TRUE, otherwise it will be AT_FALSE.

5.1.4 AT_ISREADONLY

```
int AT_IsReadOnly(AT_H Hndl, AT_WC* Feature, AT_BOOL* ReadOnly)
```

Description

This function can be used to determine whether the feature specified by the Feature parameter can be modified. On return the ReadOnly parameter will contain the value AT_FALSE or AT_TRUE. In the case that the feature can be modified the value of ReadOnly will be AT_TRUE, otherwise it will be AT_FALSE.

5.1.5 AT_IsWRITABLE

```
int AT_IsWritable(AT_H Hndl, AT_WC* Feature, AT_BOOL* Writable)
```

Description

This function can be used to determine whether the feature specified by the Feature parameter can currently be modified. On return the Writable parameter will contain the value AT_FALSE or AT_TRUE. In the case that the feature is currently writable the value of Writable will be AT_TRUE, otherwise it will be AT_FALSE. This function differs from the AT_IsReadOnly function in that a feature that is not writable may only be temporarily unavailable for modification because of the values of other features, whereas a feature that is read only is permanently un-modifiable.

5.1.6 AT_IsREADABLE

```
int AT_IsReadable(AT_H Hndl, AT_WC* Feature, AT_BOOL* Readable)
```

Description

This function can be used to determine whether the feature specified by the Feature parameter can currently be read. On return the Readable parameter will contain the value AT_FALSE or AT_TRUE. In the case that the feature is currently readable the value of Readable will be AT_TRUE, otherwise it will be AT_FALSE. A feature may become unavailable for reading based on the value of other features.

5.1.7 AT_REGISTERFEATURECALLBACK

```
int AT_RegisterFeatureCallback(AT_H Hndl, AT_WC* Feature,  
    FeatureCallback EvCallback, void* Context)
```

Description

To retrieve a notification each time the value or other properties of a feature changes you can use this function to register a callback function. The Feature that you wish to receive notifications for is passed into the function along with the function that you wish to get called. The fourth parameter is a caller defined parameter that can be used to provide contextual information when the callback is called. The callback function should have the signature shown below.

```
int AT_EXP_CONV MyFunction(AT_H Hndl, AT_WC* Feature, void* Context)
```

When called, the Feature that caused the callback is returned which allows you to use a single callback function to handle multiple features. The context parameter is the same as that used when registering the callback and is sent unmodified. As soon as this callback is registered a single callback will be made immediately to allow the callback handling code to perform any Initialisation code to set up monitoring of the feature.

5.1.8 AT_UNREGISTERFEATURECALLBACK

```
int AT_UnregisterFeatureCallback(AT_H Hndl, AT_WC* Feature,  
    FeatureCallback EvCallback, void* Context)
```

Description

This function is used to un-register a callback function previously registered using AT_RegisterFeatureCallback. The same parameters that were passed to the register function should be passed to this unregister function. Once this function is called, no more callbacks will be sent to this callback function for the specified Feature.

5.1.9 AT_INITIALISELIBRARY

```
int AT_InitialiseLibrary()
```

Description

This function is used to prepare the SDK internal structures for use and must be called before any other SDK functions have been called.

5.1.10 AT_FINALISELIBRARY

```
int AT_FinaliseLibrary()
```

Description

This function will free up any resources used by the SDK and should be called whenever the program no longer needs to use any SDK functions. AT_InitialiseLibrary may be called again later by the same process if camera control is again required.

5.1.11 AT_SETINT

```
int AT_SetInt(AT_H Hndl, AT_WC* Feature, AT_64 Value)
```

Description

This function will modify the value of the specified feature, if the feature is of integer type. The function will return an error if the feature is read only or currently not writable or if the feature is either not an integer feature or is not implemented by the camera.

5.1.12 – AT_GETINT

```
int AT_GetInt(AT_H Hndl, AT_WC* Feature, AT_64 * Value)
```

Description

This function will return the current value for the specified feature. The function will return an error if the feature is currently not readable or if the specified feature is either not an integer feature or is not implemented by the camera.

5.1.13 – AT_GETINTMAX

```
int AT_GetIntMax(AT_H Hndl, AT_WC* Feature, AT_64 * MaxValue)
```

Description

This function will return the maximum allowable value for the specified integer type feature.

5.1.14 – AT_GETINTMIN

```
int AT_GetIntMin(AT_H Hndl, AT_WC* Feature, AT_64 * MinValue)
```

Description

This function will return the minimum allowable value for the specified integer type feature.

5.1.15 – AT_SETFLOAT

```
int AT_SetFloat(AT_H Hndl, AT_WC* Feature, double Value)
```

Description

This function will modify the value of the specified feature, if the feature is of float type. The function will return an error if the feature is read only or currently not writable or if the feature is either not a float type feature or is not implemented by the camera.

5.1.16 – AT_GETFLOAT

```
int AT_GetFloat(AT_H Hndl, AT_WC* Feature, double * Value)
```

Description

This function will return the current value for the specified feature. The function will return an error if the feature is currently not readable or if the specified feature is either not a float type feature or is not implemented by the camera.

5.1.17 – AT_GETFLOATMAX

```
int AT_GetFloatMax(AT_H Hndl, AT_WC* Feature, double * MaxValue)
```

Description

This function will return the maximum allowable value for the specified float type feature.

5.1.18 – AT_GETFLOATMIN

```
int AT_GetIntMin(AT_H Hndl, AT_WC* Feature, double * MinValue)
```

Description

This function will return the minimum allowable value for the specified float type feature.

5.1.19 – AT_SETBOOL

```
int AT_SetBool(AT_H Hndl, AT_WC* Feature, AT_BOOL Value)
```

Description

This function will set the value of the specified boolean feature. A value of AT_FALSE indicates false and a value of AT_TRUE indicates true. An error will be returned if the feature is read only, currently not writable, not a boolean feature or is not implemented by the camera.

5.1.20 – AT_GETBOOL

```
int AT_GetBool(AT_H Hndl, AT_WC* Feature, AT_BOOL * Value)
```

Description

This function will return the current value of the specified boolean feature. If a value of AT_FALSE is returned then the feature is currently set to false. If a value of AT_TRUE is returned then the feature is currently set to true. An error will be returned if the feature is currently not readable, not a boolean feature or is not implemented by the camera.

5.1.21 – AT_COMMAND

```
int AT_Command(AT_H Hndl, AT_WC* Feature)
```

Description

This function will trigger the specified command feature to execute. An error will be returned if the feature is currently not writable, not a command feature or is not implemented by the camera.

5.1.22 – AT_SETSTRING

```
int AT_SetString(AT_H Hndl, AT_WC* Feature, AT_WC* Value)
```

Description

This function will set the value of the specified string feature. The string should be null terminated. An error will be returned if the feature is read only, currently not writable, not a string feature or is not implemented by the camera.

5.1.23 – AT_GETSTRING

```
int AT_GetString(AT_H Hndl, AT_WC* Feature, AT_WC* Value, int StringLength)
```

Description

This function will return the current value of the specified string feature. The length of the string in which you want the value returned must be provided in the fourth parameter and the string should include enough space for the null terminator. An error will be returned if the feature is currently not readable, not a string feature or is not implemented by the camera.

5.1.24 – AT_GETSTRINGMAXLENGTH

```
int AT_GetStringMaxLength(AT_H Hndl, AT_WC* Feature, int* MaxStringLength)
```

Description

This function will return the maximum length of the specified string feature. This value can be used to determine what size of string to allocate when retrieving the value of the feature using the AT_GetString function.

5.1.25 – AT_SETENUMINDEX

```
int AT_SetEnumIndex(AT_H Hndl, AT_WC* Feature, int Value)
```

Description

This function sets the currently selected index of the specified enumerated feature. The index is zero based and should be in the range 0 to Count-1, where Count has been retrieved using the AT_GetEnumCount function. An error will be returned if the feature is read only, currently not writable, the index is outside the allowed range, not an enumerated feature, or the feature is not implemented by the camera. In some cases an index within the range may not be allowed if its availability depends on other features values, in this case an error will be returned if this index is applied.

5.1.26 – AT_SETENUMSTRING

```
int AT_SetEnumString(AT_H Hndl, AT_WC* Feature, AT_WC* String)
```

Description

This function directly sets the current value of the specified enumerated feature. The String parameter must be one of the allowed values for the feature and must be currently available.

5.1.27 – AT_GETENUMINDEX

```
int AT_GetEnumIndex(AT_H Hndl, AT_WC* Feature, int* Value)
```

Description

This function retrieves the currently selected index of the specified enumerated feature. The function will return an error if the feature is currently not readable or if the specified feature is either not an enumerated type feature or is not implemented by the camera.

5.1.28 – AT_GETENUMCOUNT

```
int AT_GetEnumCount(AT_H Hndl, AT_WC* Feature, int* Count)
```

Description

This function returns the number of indexes that the specified enumerated feature can be set to.

5.1.29 – AT_GETENUMSTRINGBYINDEX

```
int AT_GetEnumStringByIndex(AT_H Hndl, AT_WC* Feature, int Index, AT_WC* String,  
                             int StringLength)
```

Description

This function returns the text representation of the specified enumerated feature index. The index should be in the range 0... Count-1, where Count has been retrieved using the AT_GetEnumCount function. The length of the String parameter should be passed in to the fifth parameter.

5.1.30 – AT_IsEnumIndexAvailable

```
int AT_IsEnumIndexAvailable(AT_H Hndl, AT_WC* Feature, int Index,  
                           AT_BOOL* Available)
```

Description

This function indicates whether the specified enumerated feature index can currently be selected. The availability of enumerated options may depend on the value of other features.

5.1.31 – AT_IsEnumIndexImplemented

```
int AT_IsEnumIndexImplemented(AT_H Hndl, AT_WC* Feature, int Index,  
                             AT_BOOL* Implemented)
```

Description

This function indicates whether the camera supports the specified enumerated feature index. For consistency across the camera range, some enumerated features options may appear in the list even when they are not supported, this function will let you filter out these options.

5.1.32 – AT_QueueBuffer

```
int AT_QueueBuffer(AT_H Hndl, AT_U8* Ptr, int PtrSize)
```

Description

This function configures the area of memory into which acquired images will be stored. You can call this function multiple times to set up storage for consecutive images in a series. The order in which buffers are queued is the order in which they will be used on a first in, first out (FIFO) basis. The PtrSize parameter should be equal to the size of an individual image in number of bytes. This function may be called before the acquisition starts, after the acquisition starts or a combination of the two. Any buffers queued using this function should not be modified or deallocated by the calling application until they are either returned from the AT_WaitBuffer function, or the AT_Flush function is called.

5.1.33 – AT_WaitBuffer

```
int AT_WaitBuffer(AT_H Hndl, AT_U8** Ptr, int* PtrSize, unsigned int Timeout)
```

Description

This function is used to receive notification whenever a previously queued image buffer contains data. The address of the buffer that is now available is returned in the Ptr parameter. The PtrSize parameter will return with the size of the returned image buffer. The Timeout parameter can be specified to indicate how long in milliseconds you wish to wait for the next available image. The function will put the calling thread to sleep until either an image becomes available or the Timeout elapses.

5.1.34 – AT_Flush

```
int AT_Flush(AT_H Hndl)
```

Description

This function is used to flush out any remaining buffers that have been queued using the AT_QueueBuffer function. It should always be called after the AT_Command(L"AcquisitionStop") function has been called. If this function is not called after an acquisition is complete, then the remaining buffers will be used the next time an acquisition is started, and may lead to undefined behaviour.

SECTION 6 ADDITIONAL LIBRARIES

This section describes the additional libraries that are provided as part of the SDK3 installation. These additional libraries are not required to use the SDK; they provide additional ease of use functionality.

6.1 ATUTILITY

This library provides additional general utility functions and can be used with an Embarcadero or Microsoft compatible compiler. To use this library perform the following steps (assumes you have already setup your project as described in **Section 1.6 Getting Started**.

1. Add "atutility.h" to the list of header files included in your application source file.
2. Add the appropriate library from the SDK3 installation directory to your project.
 - atutility.lib for the Embarcadero compiler
 - atutilitym.lib for the Microsoft compiler
3. Copy the "atutility.dll" from the SDK3 installation directory to the directory that the executable is going to run from.

6.1.1 – AT_INITIALISEUTILITYLIBRARY

```
int AT_InitialiseUtilityLibrary ()
```

This function is used to initialize the utility library. It must be called before any utility functions can be called. There are no parameters expected.

The following is a brief explanation of the error codes that can be returned by the function:

Error Code	Description
AT_SUCCESS (0)	The library has been initialised successfully.

6.1.2 – AT_FINALISEUTILITYLIBRARY

```
int AT_FinaliseUtilityLibrary ()
```

This function is used to close the utility library. It must be called before the user application completes to clean up internal resources. There are no parameters expected.

The following is a brief explanation of the error codes that can be returned by the function:

Error Code	Description
AT_SUCCESS (0)	The library has been closed successfully.

6.1.3 – AT_CONVERTBUFFER

```
int AT_ConvertBuffer(AT_U8* inputBuffer, AT_U8* outputBuffer, AT_64 width, AT_64 height, AT_64 stride, const AT_WC* inputPixelEncoding, const AT_WC* outputPixelEncoding)
```

This function is used to convert a buffer from one pixel encoding to another pixel encoding. For an explanation of the different encoding types see the pixel encoding feature description in the **Feature Reference Section 4.2 Feature Reference**. The converted image data will not contain any metadata therefore, if you require access to the metadata you will have to extract this from the input buffer as described in the **Metadata Section 4.5 Metadata**. The following table provides a brief description of the parameters.

Parameter	Description
inputBuffer	This is a pointer to the input buffer that you want to convert.

outputBuffer	This is a pointer to the buffer that you want the converted data to be stored in. This should be large enough to hold the converted image i.e. width x height x bytes/pixel for the output pixel encoding (Mono16 = 2, Mono32 = 4).
width	This is the width of the image stored in the input buffer in pixels. Can be determined by using the SDK3 integer feature "AOI Width".
height	This is the height of the image stored in the input buffer in pixels. Can be determined by using the SDK3 integer feature "AOI Height".
stride	This is the number of bytes/line for the image stored in the input buffer. Can be determined by using the SDK3 integer feature "AOI Stride".
inputPixelEncoding	This is the pixel encoding that was used to create the image stored in the input buffer. The valid values that can be used are: <ul style="list-style-type: none"> • Mono12 • Mono12Packed • Mono16 • Mono32
outputPixelEncoding	This is the pixel encoding that will be used to store the image in the output buffer. The valid values that can be used are: <ul style="list-style-type: none"> • Mono16 • Mono32

The following table provides a brief description of the error codes that can be returned by the function:

Error Code	Description
AT_SUCCESS (0)	The function call has been successful.
AT_ERR_NOTINITIALISED (1)	The library has not been initialised.
AT_ERR_INVALIDINPUTPIXELENCODING (1003)	The input pixel encoding is not valid.
AT_ERR_INVALIDOUTPUTPIXELENCODING (1002)	The output pixel encoding is not valid.

6.1.3.1 CONVERT BUFFER EXAMPLE

The following simple console application shows how to use the AT_ConvertBuffer function to convert from "Mono12Packed" image data to "Mono16" image data:

```
#include "atcore.h"
#include "atutility.h"

int main(int argc, char* argv[])
{
    int i_retCode;
    i_retCode = AT_InitialiseLibrary();
    if (i_retCode == AT_SUCCESS) {
        i_retCode = AT_InitialiseUtilityLibrary ();
        if (i_retCode == AT_SUCCESS) {
            AT_64 iNumberDevices = 0;
            AT_GetInt(AT_HANDLE_SYSTEM, L"Device Count", &iNumberDevices);
            if (iNumberDevices > 0) {
                AT_H Hndl;
                i_retCode = AT_Open(0, &Hndl);
            }
        }
    }
}
```

```

if (i_retCode == AT_SUCCESS) {
    AT_SetEnumeratedString(Hndl, L"Pixel Encoding", L"Mono12Packed");
    AT_SetFloat(Hndl, L"Exposure Time", 0.01);

    double temperature = 0;
    AT_SetBool(Hndl, L"SensorCooling", AT_TRUE);

    int temperatureCount = 0;
    AT_GetEnumCount(Hndl, L"TemperatureControl", &temperatureCount);
    AT_SetEnumIndex(Hndl, L"TemperatureControl", temperatureCount-1);

    int temperatureStatusIndex = 0;
    wchar_t temperatureStatus[256];
    AT_GetEnumIndex(Hndl, L"TemperatureStatus", &temperatureStatusIndex);
    AT_GetEnumStringByIndex(Hndl, L"TemperatureStatus", temperatureStatusIndex,
temperatureStatus, 256);

    while(wcsncmp(L"Stabilised",temperatureStatus) != 0) {
        AT_GetEnumIndex(Hndl, L"TemperatureStatus", &temperatureStatusIndex);
        AT_GetEnumStringByIndex(Hndl, L"TemperatureStatus", temperatureStatusIndex,
temperatureStatus, 256);
    }

    //Get the number of bytes required to store one frame
    AT_64 iImageSizeBytes;
    AT_GetInt(Hndl, L"Image Size Bytes", &iImageSizeBytes);
    int iBufferSize = static_cast<int>(iImageSizeBytes);

    //Allocate a memory buffer to store one frame
    unsigned char* UserBuffer = new unsigned char[iBufferSize];

    AT_QueueBuffer(Hndl, UserBuffer, iBufferSize);
    AT_Command(Hndl, L"Acquisition Start");

    unsigned char* Buffer;
    if (AT_WaitBuffer(Hndl, &Buffer, &iBufferSize, 10000) == AT_SUCCESS){
        //Unpack the 12 bit packed data
        AT_64 ImageHeight;
        AT_GetInt(Hndl, L"AOI Height", &ImageHeight);
        AT_64 ImageWidth;
        AT_GetInt(Hndl, L"AOI Width", &ImageWidth);
        AT_64 ImageStride;
        AT_GetInt(Hndl, L"AOI Stride", &ImageStride);
        unsigned short* unpackedBuffer =
            new unsigned short[ImageHeight*ImageWidth];
        AT_ConvertBuffer(Buffer,
            reinterpret_cast<unsigned char*>(unpackedBuffer),
            ImageWidth, ImageHeight,
            ImageStride, L"Mono12Packed", L"Mono16");

        // process unpacked image data

        delete[] unpackedBuffer;
    }
    AT_Command(Hndl, L"Acquisition Stop");
    AT_Flush(Hndl);
    delete[] UserBuffer;
}
AT_Close(Hndl);
}

```

```

    }
}
AT_FinaliseLibrary();
AT_FinaliseUtilityLibrary();
return 0;
}

```

6.1.4 – AT_CONVERTBUFFERUSINGMETADATA

```

int AT_ConvertBufferUsingMetaData(AT_U8* inputBuffer, AT_U8* outputBuffer, AT_64
imagesizebytes, const AT_WC* outputPixelEncoding)

```

This function is used to convert a buffer from the input pixel encoding to another pixel encoding. For an explanation of the different encoding types see the pixel encoding feature description in the **Feature Reference Section 4.2** Feature Reference. The converted image data will not contain any metadata therefore, if you require access to the metadata you will have to extract this from the input buffer as described in the **Metadata Section 4.5** Metadata. The following table provides a brief description of the parameters.

Parameter	Description
inputBuffer	This is a pointer to the input buffer that you want to convert.
outputBuffer	This is a pointer to the buffer that you want the converted data to be stored in. This should be large enough to hold the converted image i.e. width x height x bytes/pixel for the output pixel encoding (Mono16 = 2, Mono32 = 4).
imagesizebytes	This is the size of the image stored in the input buffer in bytes. Can be determined by using the SDK3 integer feature "ImageSizeBytes".
outputPixelEncoding	This is the pixel encoding that will be used to store the image in the output buffer. The valid values that can be used are: <ul style="list-style-type: none"> • Mono16 • Mono32

The following table provides a brief description of the error codes that can be returned by the function:

Error Code	Description
AT_SUCCESS (0)	The function call has been successful.
AT_ERR_NOTINITIALISED (1)	The library has not been initialised.
AT_ERR_INVALIDOUTPUTPIXELENCODING (1002)	The output pixel encoding is not valid.
AT_ERR_INVALIDMETADATAINFO (1004)	The input buffer does not include metadata. This may be due to the system not supporting this option or it not being activated.

6.1.5 – AT_GETWIDTHFROMMETADATA

```

int AT_GetWidthFromMetadata(AT_U8* inputBuffer, AT_64 imagesizebytes, AT_64& width)

```

This function is used to retrieve the width of the image stored in the input buffer from the metadata of the image. For more information on the metadata, see the **Metadata Section 4.5** Metadata. The following table provides a brief description of the parameters.

Parameter	Description
inputBuffer	This is a pointer to the input buffer that you want to convert.
imagesizebytes	This is the size of the image stored in the input buffer in bytes. Can be determined by using the SDK3 integer feature "ImageSizeBytes".

width	This is a reference to the variable that will contain the width of the image stored in the input buffer in pixels.
--------------	--

The following table provides a brief description of the error codes that can be returned by the function:

Error Code	Description
AT_SUCCESS (0)	The function call has been successful.
AT_ERR_METADATAANOTFOUND (1006)	The input buffer does not include the required metadata. This may be due to the system not supporting this option or it not being activated.

6.1.6 – AT_GETHEIGHTFROMMETADATA

```
int AT_GetHeightFromMetadata(AT_U8* inputBuffer, AT_64 imagesizebytes, AT_64& height)
```

This function is used to retrieve the height of the image stored in the input buffer from the metadata of the image. For more information on the metadata, see the **Metadata Section 4.5 Metadata**. The following table provides a brief description of the parameters.

Parameter	Description
inputBuffer	This is a pointer to the input buffer that you want to convert.
imagesizebytes	This is the size of the image stored in the input buffer in bytes. Can be determined by using the SDK3 integer feature "ImageSizeBytes".
height	This is a reference to the variable that will contain the height of the image stored in the input buffer in pixels.

The following table provides a brief description of the error codes that can be returned by the function:

Error Code	Description
AT_SUCCESS (0)	The function call has been successful.
AT_ERR_METADATAANOTFOUND (1006)	The input buffer does not include the required metadata. This may be due to the system not supporting this option or it not being activated.

6.1.7 – AT_GETSTRIDEFROMMETADATA

```
int AT_GetStrideFromMetadata(AT_U8* inputBuffer, AT_64 imagesizebytes, AT_64& stride)
```

This function is used to retrieve the number of bytes/line for the image stored in the input buffer from the metadata of the image. For more information on the metadata, see the **Metadata Section 4.5 Metadata**. The following table provides a brief description of the parameters.

Parameter	Description
inputBuffer	This is a pointer to the input buffer that you want to convert.
imagesizebytes	This is the size of the image stored in the input buffer in bytes. Can be determined by using the SDK3 integer feature "ImageSizeBytes".
stride	This is a reference to the variable that will contain the number of bytes/line for the image stored in the input buffer.

The following table provides a brief description of the error codes that can be returned by the function:

Error Code	Description
AT_SUCCESS (0)	The function call has been successful.
AT_ERR_METADATANOTFOUND (1006)	The input buffer does not include the required metadata. This may be due to the system not supporting this option or it not being activated.

6.1.8 – AT_GETPIXELENCODINGFROMMETADATA

```
int AT_GetPixelEncodingFromMetadata(AT_U8* inputBuffer, AT_64 imagesizebytes, AT_WC* pixelEncoding, AT_U8 pixelEncodingSize)
```

This function is used to retrieve the pixel encoding of the image stored in the input buffer from the metadata of the image. For an explanation of the different encoding types see the pixel encoding feature description in the **Feature Reference Section 4.2** Feature Reference. For more information on the metadata, see the **Metadata Section 4.5** Metadata. The following table provides a brief description of the parameters.

Parameter	Description
inputBuffer	This is a pointer to the input buffer that you want to convert.
imagesizebytes	This is the size of the image stored in the input buffer in bytes. Can be determined by using the SDK3 integer feature "ImageSizeBytes".
pixelEncoding	This is the pixel encoding that was used to create the image stored in the input buffer. The valid values that can be used are: <ul style="list-style-type: none"> • Mono12 • Mono12Packed • Mono16 • Mono32
pixelEncodingSize	This is the number of characters that the pixelEncoding variable can contain.

The following table provides a brief description of the error codes that can be returned by the function:

Error Code	Description
AT_SUCCESS (0)	The function call has been successful.
AT_ERR_METADATANOTFOUND (1006)	The input buffer does not include the required metadata. This may be due to the system not supporting this option or it not being activated.

6.1.9 – AT_GETTIMESTAMPFROMMETADATA

```
int AT_GetTimeStampFromMetadata(AT_U8* inputBuffer, AT_64 imagesizebytes, AT_64& timeStamp)
```

This function is used to retrieve the timestamp of the image stored in the input buffer from the metadata of the image. For more information on the metadata, see the **Metadata Section 4.5** Metadata. The following table provides a brief description of the parameters.

Parameter	Description
inputBuffer	This is a pointer to the input buffer that you want to convert.
imagesizebytes	This is the size of the image stored in the input buffer in bytes. Can be determined

	by using the SDK3 integer feature "ImageSizeBytes".
timeStamp	This is a reference to the variable that will contain the timestamp of the image stored in the input buffer in pixels.

The following table provides a brief description of the error codes that can be returned by the function:

Error Code	Description
AT_SUCCESS (0)	The function call has been successful.
AT_ERR_METADATANOTFOUND (1006)	The input buffer does not include the required metadata. This may be due to the system not supporting this option or it not being activated.

APPENDIX A sCMOS Feature Quick Reference

Feature	Type	Available Options
AccumulateCount	Integer	Na
AcquisitionStart	Command	Na
AcquisitionStop	Command	Na
AlternatingReadoutDirection	Boolean	Na
AOIBinning	Enumerated	1x1, 2x2, 3x3, 4x4, 8x8
AOIBin	Integer	Na
AOIHeight	Integer	Na
AOILayout	Enumerated	Image, Multitrack
AOILeft	Integer	Na
AOIStride	Integer	Na
AOITop	Integer	Na
AOIVBin	Integer	Na
AOIWidth	Integer	Na
AuxiliaryOutSource	Enumerated	FireRow1, FireRowN, FireAll, FireAny
AuxOutSourceTwo	Enumerated	ExternalShutterControl, FrameClock, RowClock, ExposedRowClock
Baseline	Integer	Na
BitDepth	Enumerated	11 Bit or 12Bit, 16 Bit
BufferOverflowEvent	Integer	Na
BytesPerPixel	Floating Point	Na
CameraAcquiring	Boolean	Na
CameraDump	Command	Na
CameraModel	String	Na
CameraName	String	Na
CameraPresent	Boolean	Na
ControllerID	String	Na
FrameCount	Integer	Na
CycleMode	Enumerated	Fixed, Continuous
DeviceCount	Integer	Na
ElectronicShutteringMode	Enumerated	Rolling, Global
EventEnable	Boolean	Na
EventsMissedEvent	Integer	Na
EventSelector	Enumerated	ExposureEndEvent, ExposureStartEvent, RowNExposureEndEvent, RowNExposureStartEvent, EventsMissedEvent, BufferOverflowEvent
ExposedPixelHeight	Integer	Na
ExposureTime	Floating Point	Na
ExposureEndEvent	Integer	Na
ExposureStartEvent	Integer	Na
External Trigger Delay	Floating Point	Na
FanSpeed	Enumerated	Off, Low, On
FastAOIFrameRateEnable	Boolean	Na
FirmwareVersion	String	Na
FrameRate	Floating Point	Na
FullAOIControl	Boolean	Na
ImageSizeBytes	Integer	Na
InterfaceType	String	Na
IOInvert	Boolean	Na
IOSelector	Enumerated	Fire 1, Fire N, Aux Out 1, Arm, External Trigger
LineScanSpeed	Floating Point	Na
LUTIndex	Integer	Na
LUTValue	Integer	Na
MaxInterfaceTransferRate	Floating Point	Na
MetadataEnable	Boolean	Na
MetadataTimestamp	Boolean	Na
MetadataFrame	Boolean	Na
MultitrackBinned	Boolean	Na
MultitrackCount	Integer	Na

Feature	Type	Available Options
MultitrackEnd	Integer	Na
MultitrackSelector	Integer	Na
MultitrackStart	Integer	Na
Overlap	Boolean	Na
PixelEncoding	Enumerated	Mono12, Mono12Packed, Mono16, Mono32
PixelHeight	Floating Point	Na
PixelReadoutRate	Enumerated	280 MHz, 200 MHz, 100 MHz
PixelWidth	Floating Point	Na
PreAmpGainControl	Enumerated	Gain 1 (11 bit), Gain 2 (11 bit), Gain 3 (11 bit), Gain 4 (11 bit), Gain 1 Gain 3 (16 bit), Gain 1 Gain 4 (16 bit), Gain 2 Gain 3 (16 bit), Gain 2 Gain 4 (16 bit)
ReadoutTime	Floating Point	Na
RollingShutterGlobalClear	Boolean	Na
RowNExposureEndEvent	Integer	Na
RowNExposureStartEvent	Integer	Na
RowReadTime	Floating Point	Na
ScanSpeedControlEnable	Boolean	Na
SensorCooling	Boolean	Na
SensorHeight	Integer	Na
SensorReadoutMode	Enumerated	Bottom Up Sequential, Bottom Up Simultaneous, Centre Out Simultaneous, Outside In Simultaneous, Top Down Sequential, Top Down Simultaneous
SensorTemperature	Floating Point	Na
SensorWidth	Integer	Na
SerialNumber	String	Na
ShutterOutputMode	Enumerated	Open, Closed
SimplePreAmpGainControl	Enumerated	11-bit (high well capacity) or 12-bit (high well capacity), 11-bit (low noise) or 12-bit (low noise), 16-bit (low noise & high well capacity)
Shutter Transfer Time	Floating Point	Na
SoftwareTrigger	Command	Na
StaticBlemishCorrection	Boolean	Na
SpuriousNoiseFilter	Boolean	Na
TargetSensorTemperature	Floating Point	Na
TemperatureControl	Enumerated	Na
TemperatureStatus	Enumerated	Cooler Off, Stabilised, Cooling, Drift, Not Stabilised, Fault
TimestampClock	Integer	Na
TimestampClockFrequency	Integer	Na
TimestampClockReset	Command	Na
TriggerMode	Enumerated	Internal, Software, External, External Start, External Exposure
VerticallyCentreAOI	Boolean	Na

APPENDIX B Apogee Feature Quick Reference

Feature	Type	Available Options
AcquiredCount	Integer	Na
AOIBin	Integer	Na
AOIHeight	Integer	Na
AOILayout	Enumerated	Image, Kinetics, TDI
AOILeft	Integer	Na
AOITop	Integer	Na
AOIVBin	Integer	Na
AOIWidth	Integer	Na
BackoffTemperatureOffset	Floating Point	Na
BitDepth	Enumerated	12 Bit(Not AltaF/Aspen/Ascent), 16 Bit For AltaU/E use PixelReadoutRate "Normal" - > "16 bit", "Fast" -> "12-bit"
CameraName	String	Na
CameraFamily	String	Na
CameraMemory	Integer	Na
ColourFilter	Enumerated	None, Blue, TrueSense
CoolerPower	Double	Na
DDR2Type	String	Na
DisableShutter	Boolean	Na
DriverVersion	String	Na
ExternalIOReadout	Boolean	Na
FanSpeed	Enumerated	Off, Low, Medium, High
FirmwareVersion	String	Na
ForceShutterOpen	String	Na
FrameCount	Integer	Na
FrameRate	Floating Point	Na
FrameInterval	Boolean	Na
FrameIntervalTiming	Boolean	Na
HeatSinkTemperature	Floating Point	Na
InputVoltage	Floating Point	Na
InterfaceType	String	Na
IOControl	Enumerated	Default, User
IODirection	Boolean/Enumerated	Input/0, Output/1
IOState	Boolean	Na
IRPreFlashEnable	Boolean	Na
KeepCleanEnable	Boolean	Na
KeepCleanPostExposureEnable	Boolean	Na
MicrocodeVersion	String	Na
Overlap	Boolean	Na
PixelHeight	Floating Point	Na
PixelReadoutRate	Enumerated	Normal, Fast (Not Available AltaE)
PortSelector	Integer	Na
PreAmpGainValue	Integer	Na
PreAmpOffsetValue	Integer	Na
SensorCooling	Boolean	Na
SensorHeight	Integer	Na
SensorModel	String	Na
SensorTemperature	Floating Point	Na
SensorType	Enumerated	CCD, CMOS
SensorWidth	Integer	Na
ShutterMode	Enumerated	Open, Closed, Auto
ShutterStrobePeriod	Double	Na
ShutterStrobePosition	Double	Na
ShutterAmpControl	Boolean	Na
ShutterState	Boolean	Na
TemperatureStatus	Enumerated	Backoff
TransmitFrames	Boolean	Na
TriggerMode	Enumerated	Internal, External, External Start, External Exposure

Feature	Type	Available Options
UsbProductId	Integer	Na
UsbDeviceId	Integer	Na

APPENDIX C Function Quick Reference

```
int AT_InitialiseLibrary();
int AT_FinaliseLibrary();

int AT_Open(int DeviceIndex, AT_H* Handle);
int AT_Close(AT_H Hndl);

typedef int (*FeatureCallback)(AT_H Hndl, AT_WC* Feature, void* Context);
int AT_RegisterFeatureCallback(AT_H Hndl, AT_WC* Feature, FeatureCallback EvCallback,
void* Context);
int AT_UnregisterFeatureCallback(AT_H Hndl, AT_WC* Feature, FeatureCallback EvCallback,
void* Context);

int AT_IsImplemented(AT_H Hndl, AT_WC* Feature, AT_BOOL* Implemented);
int AT_IsReadOnly(AT_H Hndl, AT_WC* Feature, AT_BOOL* ReadOnly);
int AT_IsReadable(AT_H Hndl, AT_WC* Feature, AT_BOOL* Readable);
int AT_IsWritable(AT_H Hndl, AT_WC* Feature, AT_BOOL* Writable);

int AT_SetInt(AT_H Hndl, AT_WC* Feature, AT_64 Value);
int AT_GetInt(AT_H Hndl, AT_WC* Feature, AT_64 * Value);
int AT_GetIntMax(AT_H Hndl, AT_WC* Feature, AT_64 * MaxValue);
int AT_GetIntMin(AT_H Hndl, AT_WC* Feature, AT_64 * MinValue);

int AT_SetFloat(AT_H Hndl, AT_WC* Feature, double Value);
int AT_GetFloat(AT_H Hndl, AT_WC* Feature, double * Value);
int AT_GetFloatMax(AT_H Hndl, AT_WC* Feature, double * MaxValue);
int AT_GetFloatMin(AT_H Hndl, AT_WC* Feature, double * MinValue);

int AT_SetBool(AT_H Hndl, AT_WC* Feature, AT_BOOL Value);
int AT_GetBool(AT_H Hndl, AT_WC* Feature, AT_BOOL* Value);

int AT_SetEnumIndex(AT_H Hndl, AT_WC* Feature, int Value);
int AT_SetEnumString(AT_H Hndl, AT_WC* Feature, AT_WC* String);
int AT_GetEnumIndex(AT_H Hndl, AT_WC* Feature, int* Value);
int AT_GetEnumCount(AT_H Hndl, AT_WC* Feature, int* Count);
int AT_IsEnumIndexAvailable(AT_H Hndl, AT_WC* Feature, int Index, AT_BOOL* Available);
int AT_IsEnumIndexImplemented(AT_H Hndl, AT_WC* Feature, int Index, AT_BOOL*
Implemented);
int AT_GetEnumStringByIndex(AT_H Hndl, AT_WC* Feature, int Index, AT_WC* String, int
StringLength);

int AT_Command(AT_H Hndl, AT_WC* Feature);

int AT_SetString(AT_H Hndl, AT_WC* Feature, AT_WC* Value);
int AT_GetString(AT_H Hndl, AT_WC* Feature, AT_WC* Value, int StringLength);
int AT_GetStringMaxLength(AT_H Hndl, AT_WC* Feature, int* MaxStringLength);

int AT_QueueBuffer(AT_H Hndl, AT_U8 * Ptr, int PtrSize);
int AT_WaitBuffer(AT_H Hndl, AT_U8 ** Ptr, int* PtrSize, unsigned int Timeout);
int AT_Flush(AT_H Hndl);
```

APPENDIX D Code Listing for Tutorial

```
#include "atcore.h"

int main(int argc, char* argv[])
{
    int i_returnCode;
    AT_H Hndl;
    int i_cameraIndex = 2;
    i_returnCode = AT_InitialiseLibrary( );

    if (i_returnCode == AT_SUCCESS) {
        i_returnCode = AT_Open ( i_cameraIndex, &Hndl );
        AT_WC ExpFeatureName[] = L"Exposure Time";
        double d_newExposure = 0.02;
        i_returnCode = AT_SetFloat ( Hndl,  ExpFeatureName, d_newExposure);
        if (i_returnCode == AT_SUCCESS) {
            //it has been set
            double d_actualExposure;
            i_returnCode = AT_GetFloat ( Hndl,  ExpFeatureName, &d_actualExposure);
            if (i_returnCode == AT_SUCCESS) {
                //the actual exposure being used is d_actualExposure

                AT_64 ImageSizeBytes;
                AT_GetInt( Hndl, L"Image Size Bytes", &ImageSizeBytes);
                //cast to prevent warnings
                int i_imageSize = static_cast<int>(ImageSizeBytes);

                unsigned char* uc_Buffer = NULL;
                unsigned char* gblp_Buffer = new unsigned char[i_imageSize+8];

                unsigned char* pucAlignedBuffer = reinterpret_cast<unsigned char*>(
                    reinterpret_cast<unsigned long>( gblp_Buffer ) + 7 ) & ~0x7);

                i_returnCode = AT_QueueBuffer(Hndl,  pucAlignedBuffer, i_imageSize);
                if (i_returnCode == AT_SUCCESS) {
                    i_returnCode = AT_Command(Hndl, L"Acquisition Start");
                    if (i_returnCode == AT_SUCCESS) {
                        unsigned char* pBuf;
                        int BufSize;
                        i_returnCode = AT_WaitBuffer(Hndl, &pBuf, &BufSize, 10000);
                        if ( i_returnCode == AT_SUCCESS) {
                            //successfully got image
                            if (pBuf == pucAlignedBuffer) {
                                //check pixel encoding to confirm format of data stream and process
                            }
                            else {
                                //Error buffer pointer incorrect from AT_WaitBuffer
                            }
                        }
                    }
                    else {
                        //error with AT_WaitBuffer, analyse i_returnCode
                    }
                }

                AT_Command(Hndl, L"Acquisition Stop");
                AT_Flush(Hndl );
            }
        }
    }
}
```

```
    }

    i_returnCode = AT_Close ( Hndl );
    if (i_returnCode != AT_SUCCESS) {
        // error closing handle
    }
}

i_returnCode = AT_FinaliseLibrary( );
if (i_returnCode != AT_SUCCESS) {
    //Error FinaliseLibrary
}
return 0;
}
```

APPENDIX E Conversion between char* and AT_WC

The following code shows an example of how to convert a char* null terminated string to the equivalent wide character string.

```
#include "stdlib.h"
char szStr[512];
AT_WC wcszStr[512];
mbstowcs(wcszStr, szStr, 512);
```

and from wide character string to char*

```
#include "stdlib.h"
char szStr[512];
AT_WC wcszStr[512];
wcstombs(szStr, wcszStr, 512);
```